

LCMS Chiral Separation of Epoxiconazole by Lux 3u Cellulose-2

Column: Lux® 3 µm Cellulose-2, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4456-B0

Elution Type: Isocratic

Eluent A: 5mM Ammonium Hydrogencarbonate

Eluent B: MeOH

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	10	90

Profile:


Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)



Products used in this application:

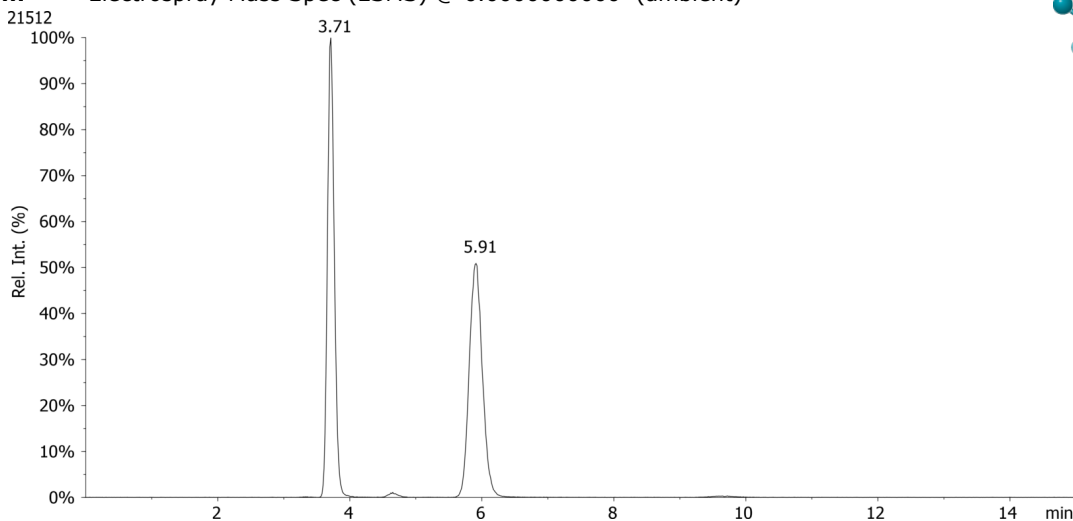

ANALYTES:

1 Epoxiconazole

Retention Time: 3.71 min

2 Epoxiconazole

Retention Time: 5.91 min

