

LCMS Chiral Separation of Triticonazole by Lux 3u Cellulose-1

Column: Lux® 3 µm Cellulose-1, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4458-B0

Elution Type: Isocratic

Eluent A: 5mM Ammonium Hydrogencarbonate

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	40	60	


Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)



Products used in this application:

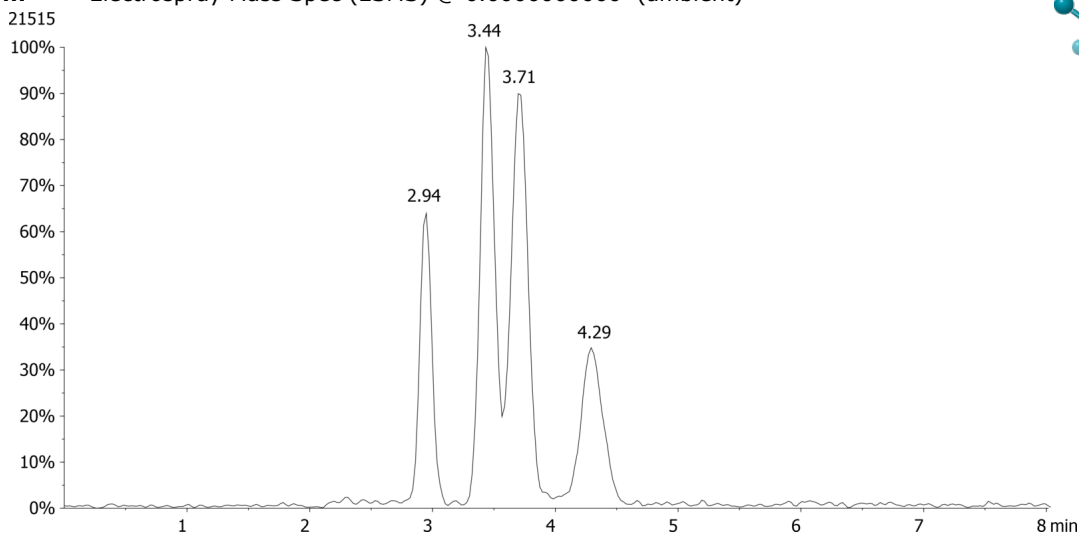

ANALYTES:

- 1 Triticonazole
Retention Time: 2.7 min
- 2 Triticonazole
Retention Time: 3.8 min
- 3 Triticonazole
Retention Time: 4.13 min
- 4 Triticonazole
Retention Time: 5.46 min

