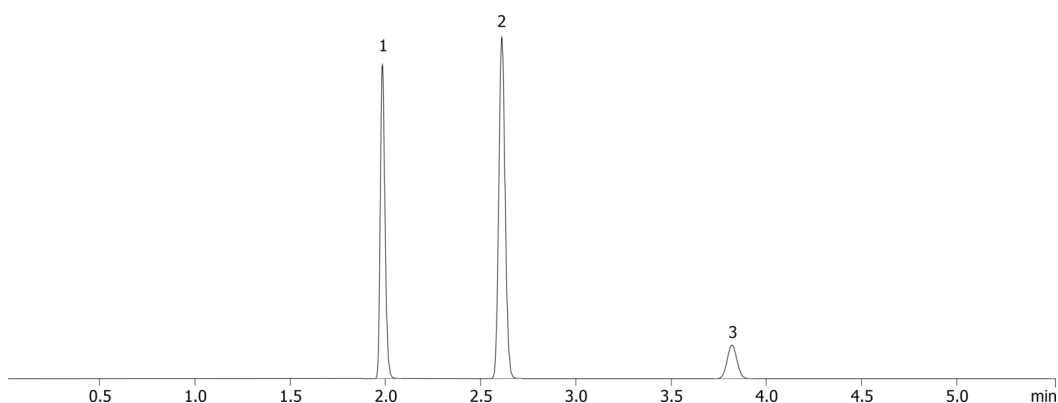


Underivatized Formaldehyde Analysis on ZB-1**Column:** Zebtron[™] ZB-1, GC Cap. Column 30 m x 0.32 mm x 3.00 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 3 µm**Order No:** 7HM-G001-36**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.5 mL/min**Injection:** Split 5.000000000:1 0 µL @ 140°C**Detection:** Refractive Index (250°C)

Analyst Note: Sample injected using a Tekmar Versa Headspace Autosampler.
Headspace oven @ 60 °C for 10 minutes; transfer line = 65 °C; Vial pressurization = 15psi.
GC analysis on Shimadzu GC-2010 plus
Analytes were 0.5% in water. Nine different samples were prepared and analyzed with 2.0% RSD in peak area for formaldehyde.

**Products used in this application:**

21516

**ANALYTES:**

- 1 Formaldehyde
- 2 Methanol
- 3 Ethanol

