

Reduced Human Ig-G1 on Aeris: LC-UV Separation

Column: Aeris™ 3.6 µm WIDEPORE XB-C8 200 Å, LC Column 100 x 2.1 mm, Ea

Dimensions: 100 x 2.1 mm ID

Order No: 00D-4481-AN

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: IPA/AcCN/0.1% TFA in water 80:10:10

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	5	60	40
	3	15	28	72
	4	19	24	76
	5	19.01	3	97
	6	20	3	97

Flow Rate: 0.4 mL/min


Col. Temp.: 70 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

21521



Products used in this application:

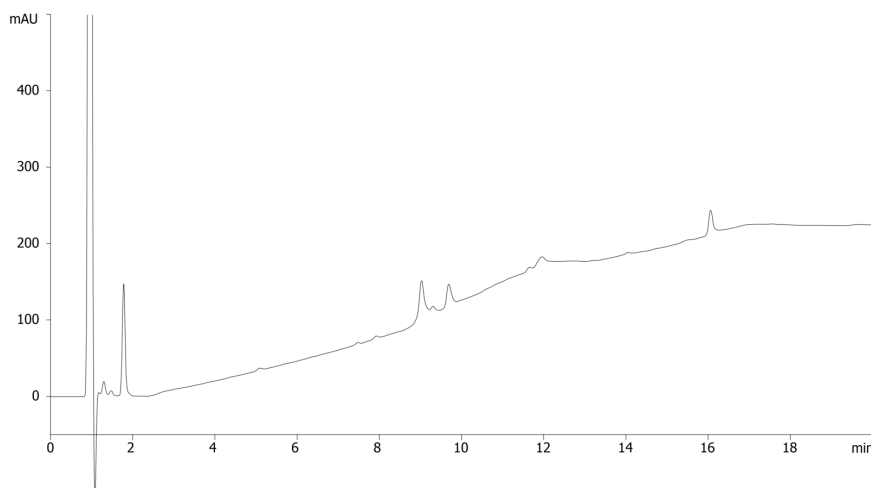

ANALYTES:

- 1 Reduced Heavy Chain Human IgG1

Retention Time: 8.9 min

- 2 Light Chain IgG1

Retention Time: 9.7 min

