


FMOC-Glu(OtBu)-OH on Lux 5 μ m Amylose-2 in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / EtOH (60:40) Formic acid (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)**ANALYTES:****1** FMOC-Glu(OtBu)-OH

Retention Time: 6.898 min

2 FMOC-Glu(OtBu)-OH

Retention Time: 12.836 min

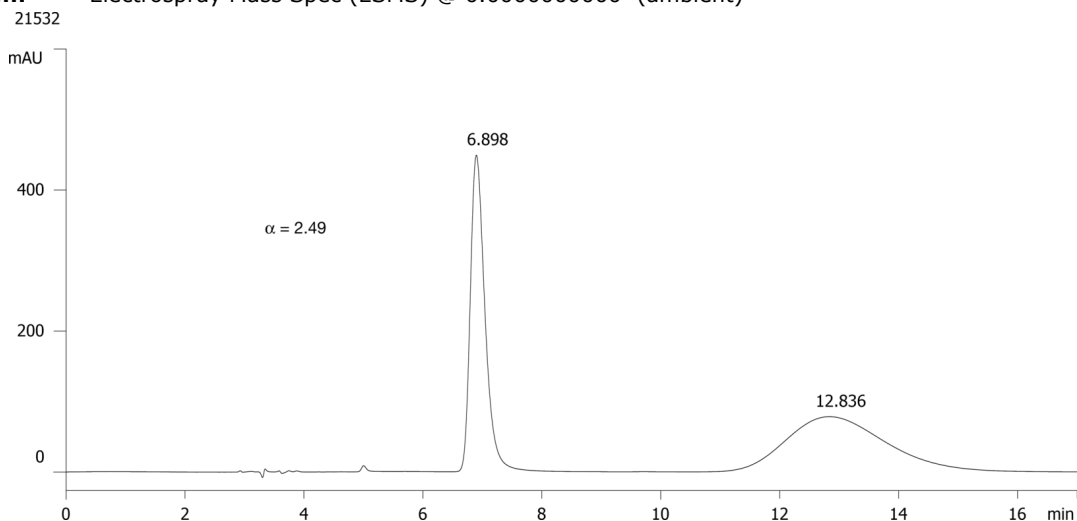