

## FMOC-Tyr(tBu)-OH on Lux 5µm Cellulose-1 in NP


**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4459-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane / IPA (60:40) Formic acid (0.1%)  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)



Products used in this application:

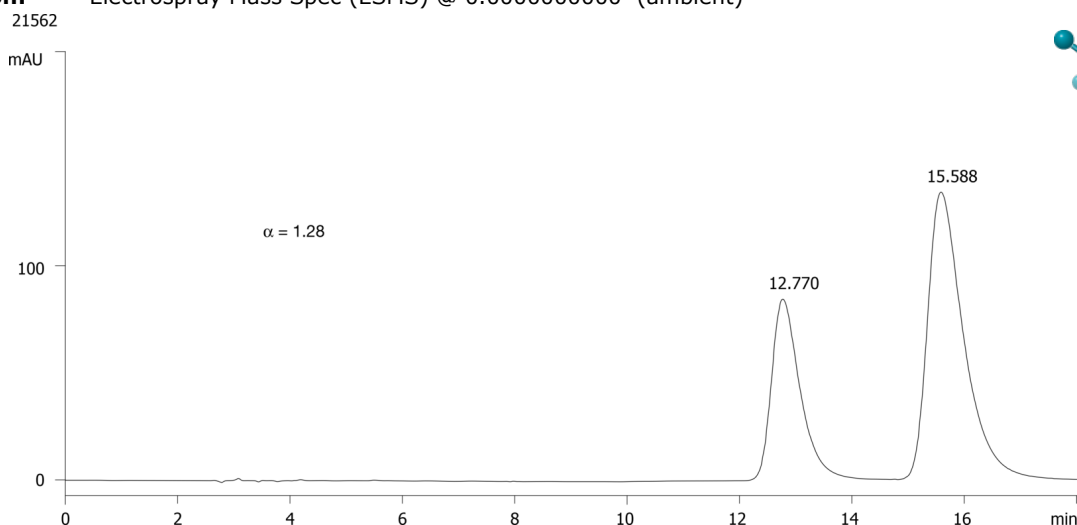

### ANALYTES:

- 1** FMOC-Tyr(tBu)-OH  
Retention Time: 12.77 min
- 2** FMOC-Tyr(tBu)-OH  
Retention Time: 15.588 min

