

## FMOC-Gln(Trt)-OH on Lux 5 $\mu$ m Cellulose-2 in RP


**Column:** Lux<sup>®</sup> 5  $\mu$ m Cellulose-2, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4457-E0  
**Elution Type:** Isocratic  
**Eluent A:** ACN / 0.1% Formic Acid (70:30) with Formic acid (0.1%) in buffer  
**Gradient**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)



Products used in this application:

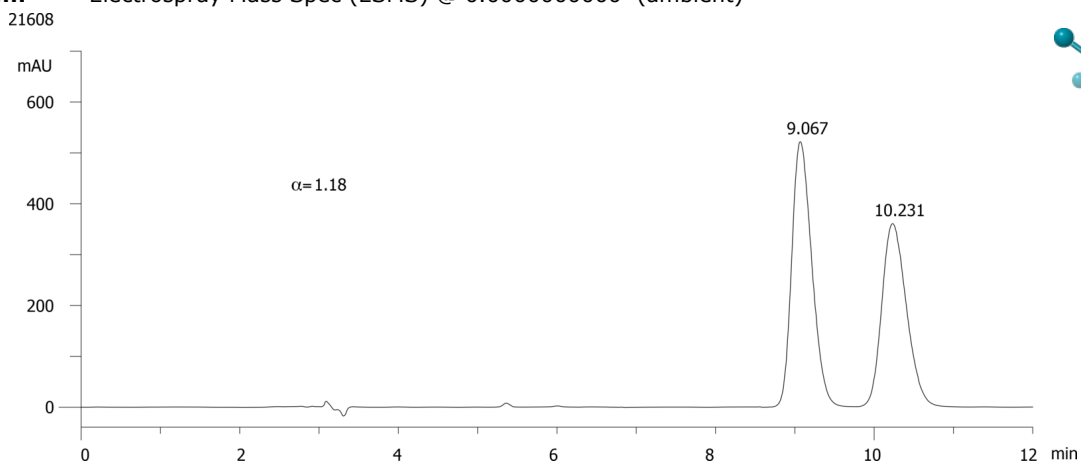

### ANALYTES:

- 1** FMOC-Gln(Trt)-OH  
Retention Time: 9.067 min
- 2** FMOC-Gln(Trt)-OH  
Retention Time: 10.231 min

