


FMOC-Glu(OtBu)-OH on Lux 5 μ m Amylose-2 in RP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** ACN / 0.1% Formic Acid (70:30) with Formic acid (0.1%) in buffer**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)**ANALYTES:****1** FMOC-Glu(OtBu)-OH

Retention Time: 4.979 min

2 FMOC-Glu(OtBu)-OH

Retention Time: 6.985 min

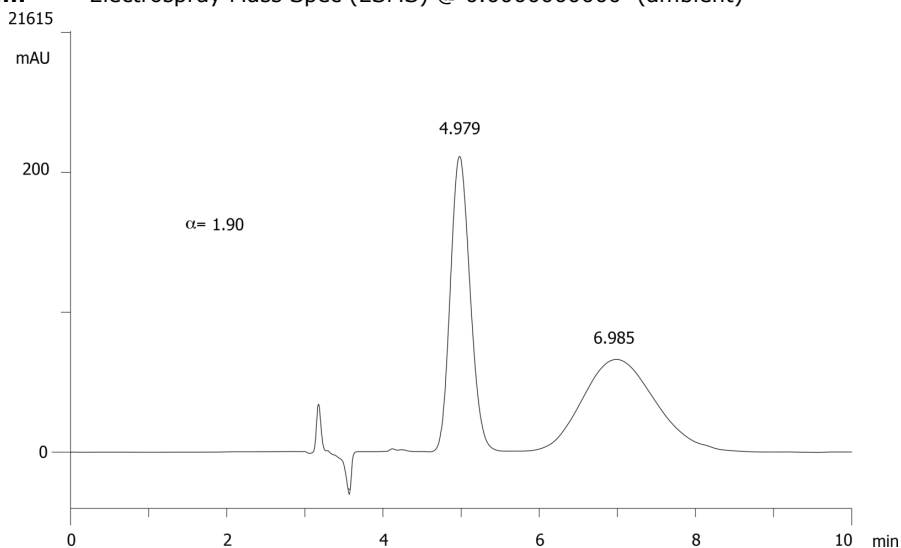