


FMOC-Val-OH on Lux 5µm Amylose-2 in RP

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: ACN / 0.1% Formic Acid (50:50) with Formic acid (0.1%) in buffer
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

21620

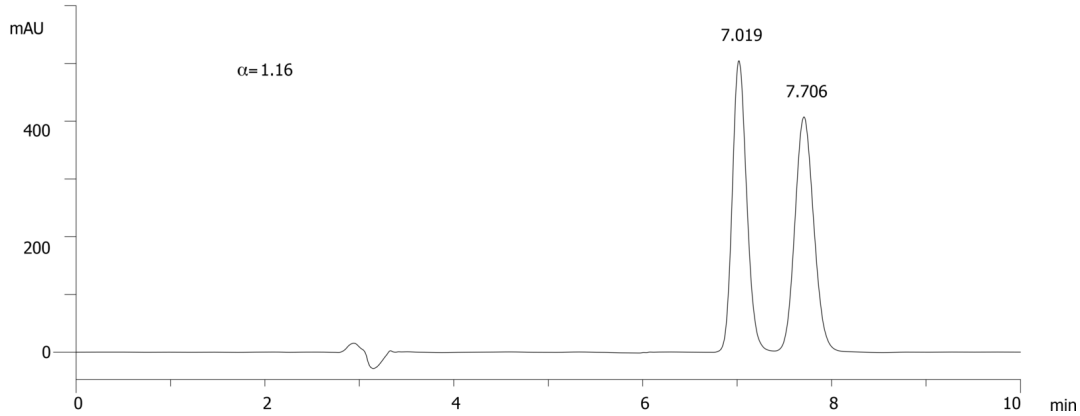

ANALYTES:

- 1 FMOC-Val-OH
Retention Time: 7.019 min
- 2 FMOC-Val-OH
Retention Time: 7.706 min

