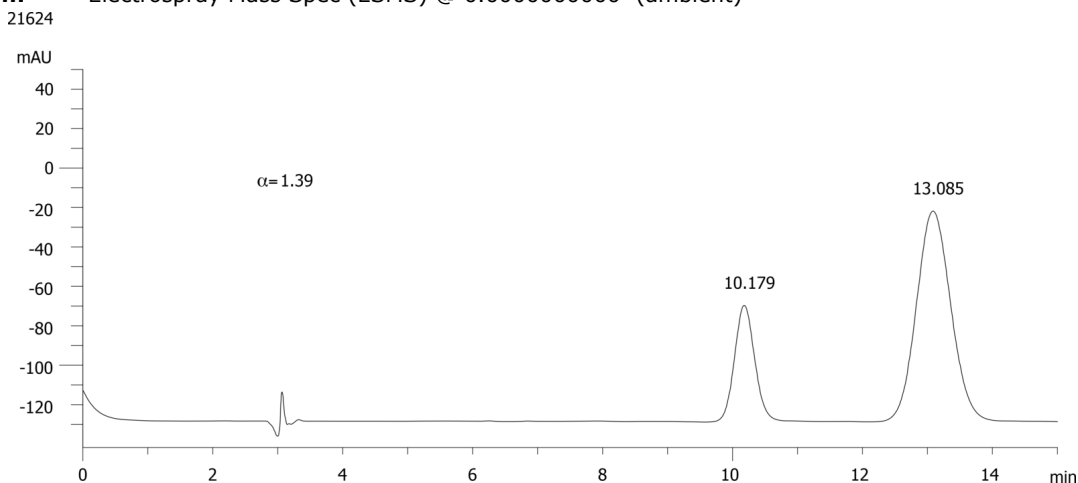


FMOC-Pro-OH on Lux 5 μ m Amylose-2 in RP

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: ACN / 0.1% TFA (40:60) with TFA (0.1%) in buffer
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

**ANALYTES:**

- 1 FMOC-Pro-OH
Retention Time: 10.179 min
- 2 FMOC-Pro-OH
Retention Time: 13.085 min

