


Diclofop methyl on Lux 5 μ m Amylose-2 in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / IPA (90:10) DEA (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)
21677**ANALYTES:****1** Diclofop methyl

Retention Time: 4.996 min

2 Diclofop methyl

Retention Time: 5.294 min

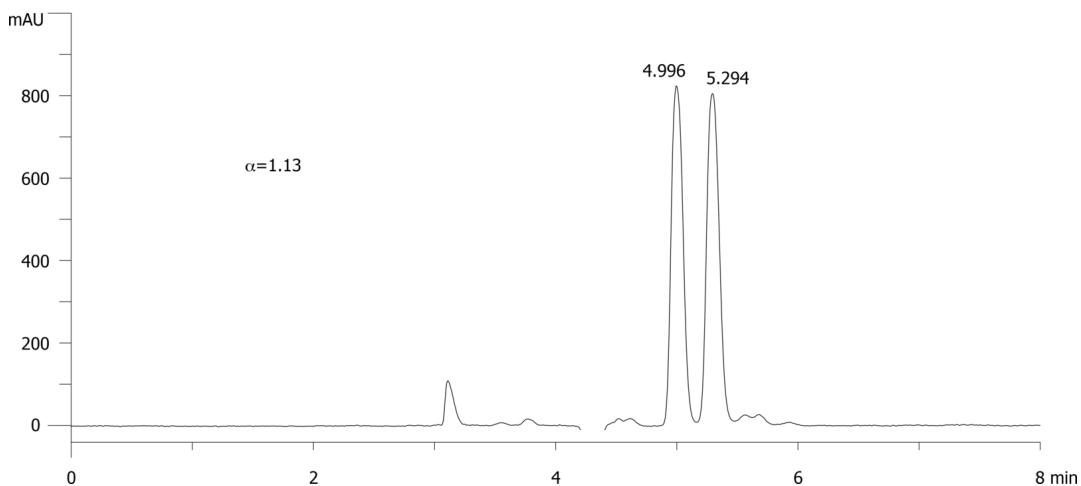