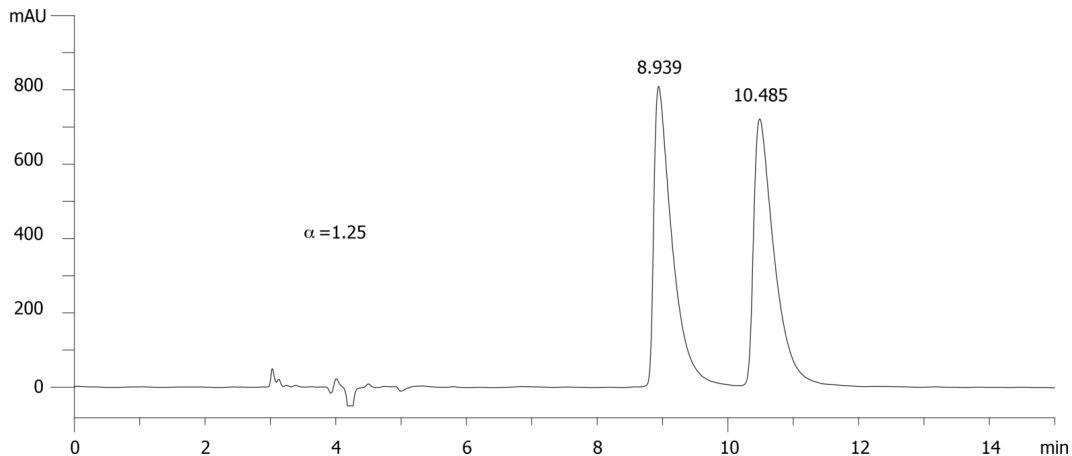


Diclofop on Lux 5 μ m Amylose-2 in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / IPA (80:20) Formic acid (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)

21687

**ANALYTES:****1** Diclofop

Retention Time: 8.939 min

2 Diclofop

Retention Time: 10.485 min

