


**Epoxiconazole on Lux 5 $\mu$ m Amylose-2 in NP****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / IPA (80:20) DEA (0.1%)**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)**ANALYTES:****1** Epoxiconazole

Retention Time: 13.667 min

**2** Epoxiconazole

Retention Time: 16.418 min

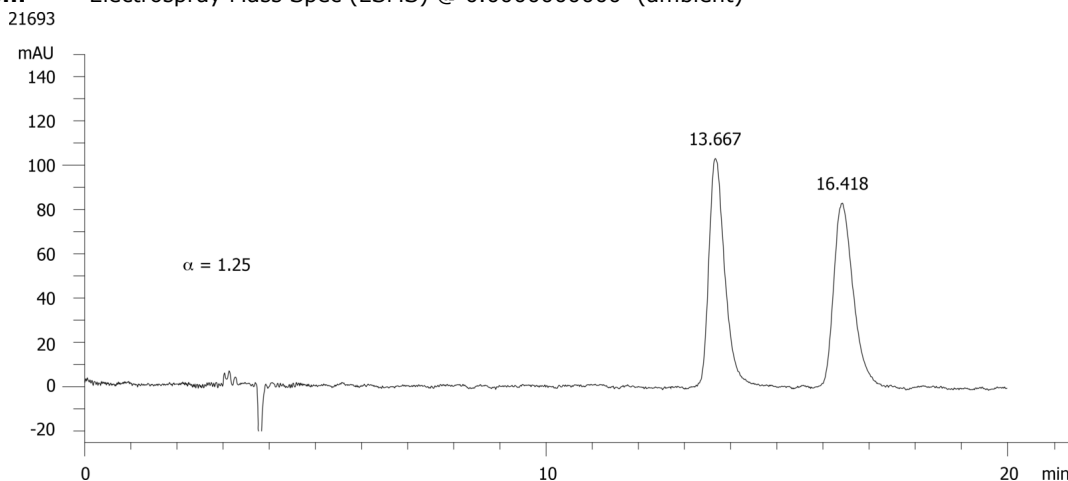