

## Fluazifop butyl on Lux 5µm Cellulose-3 in NP

**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4493-E0

**Elution Type:** Isocratic

**Eluent A:** Hexane / EtOH (95:5) DEA (0.1%)

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |


**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)



Products used in this application:

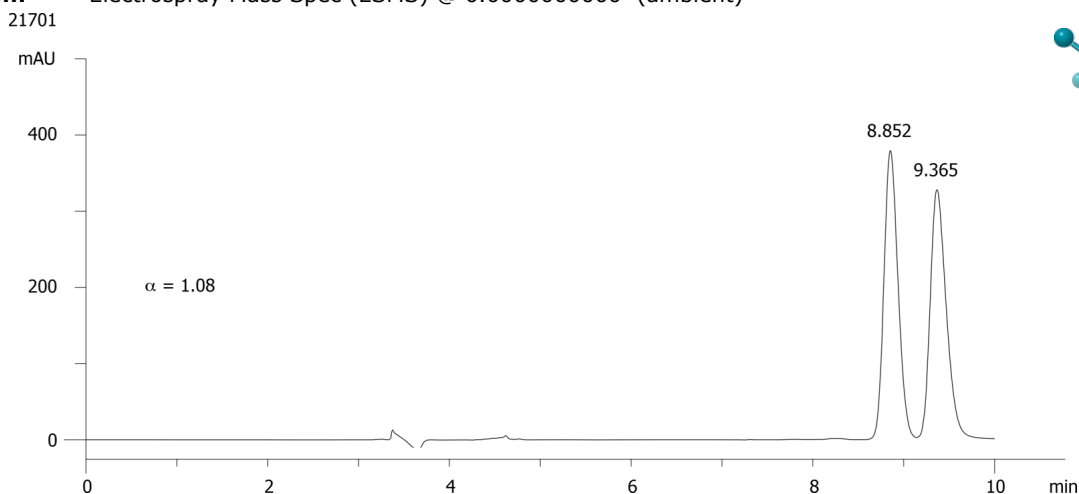

### ANALYTES:

- 1 Fluazifop butyl  
Retention Time: 8.852 min
- 2 Fluazifop butyl  
Retention Time: 9.365 min

