


Haloxyfop on Lux 5µm Amylose-2 in NP

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: Hexane / IPA (80:20) Formic acid (0.1%)
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

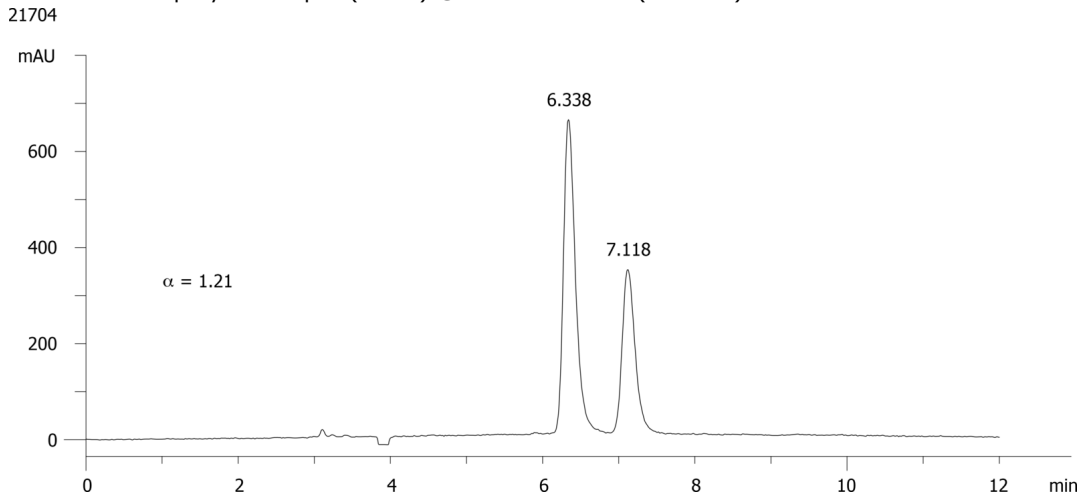

ANALYTES:

- 1 Haloxyfop
Retention Time: 6.338 min
- 2 Haloxyfop
Retention Time: 7.118 min

