


**Benalaxy on Lux 5µm Amylose-2 in RP****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** ACN / 20mM NH<sub>4</sub>HCO<sub>3</sub>+0.1%DEA (60:40) with DEA (0.1%) in buffer**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)  
App ID 21748**ANALYTES:****1** Benalaxy

Retention Time: 6.996 min

**2** Benalaxy

Retention Time: 8.087 min

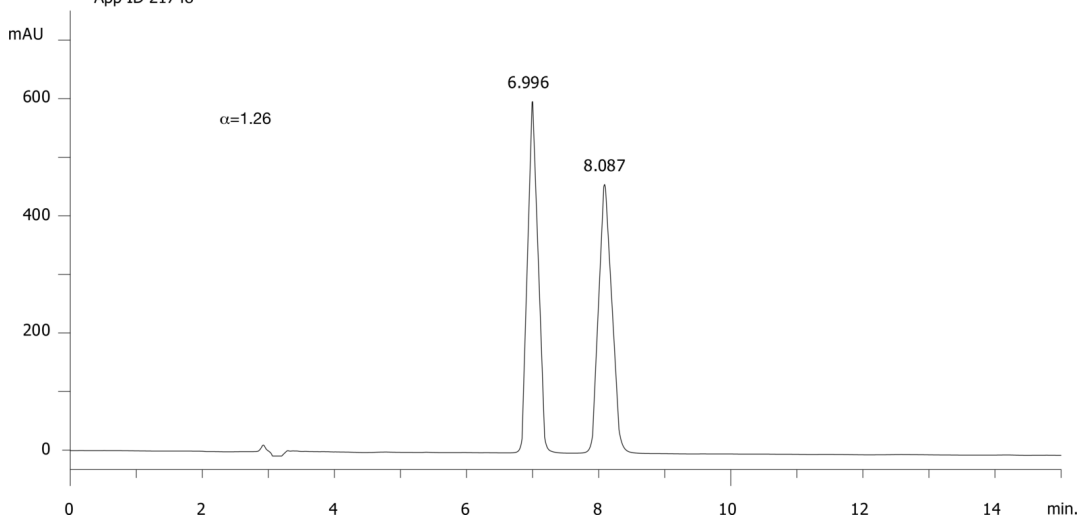