

## Diniconazole on Lux 5µm Cellulose-4 in RP


**Column:** Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4491-E0  
**Elution Type:** Isocratic  
**Eluent A:** ACN / 20mM NH<sub>4</sub>HCO<sub>3</sub>+0.1%DEA (60:40) with DEA (0.1%) in buffer  
**Gradient**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)  
 App ID 21766



Products used in this application:

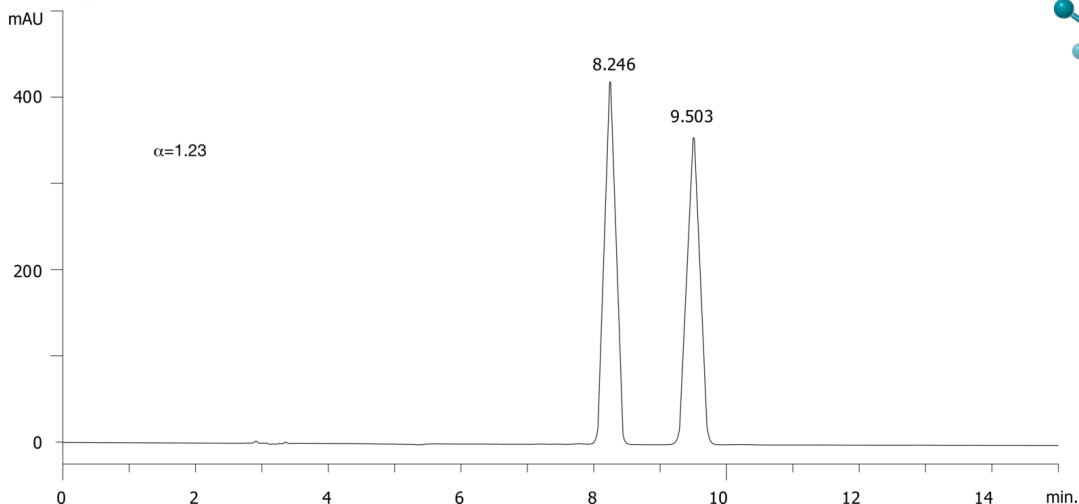

### ANALYTES:

- 1 Diniconazole  
Retention Time: 8.246 min
- 2 Diniconazole  
Retention Time: 9.503 min

