

Capsaicin and Dihydrocapsaicin Product Ion Spectra on LC/MS/MS on Kinetex 5u C18 100x2.1

Column: Kinetex® 5µm C18 100 Å, LC Column 100 x 2.1 mm, Ea

Dimensions: 100 x 2.1 mm ID

Order No: 00D-4601-AN

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in Methanol


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	7	10	90
	3	9	10	90
	4	9.1	70	30
	5	12	70	30

Flow Rate: 0.6 mL/min

Col. Temp.: 25 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (600 °C)

21863

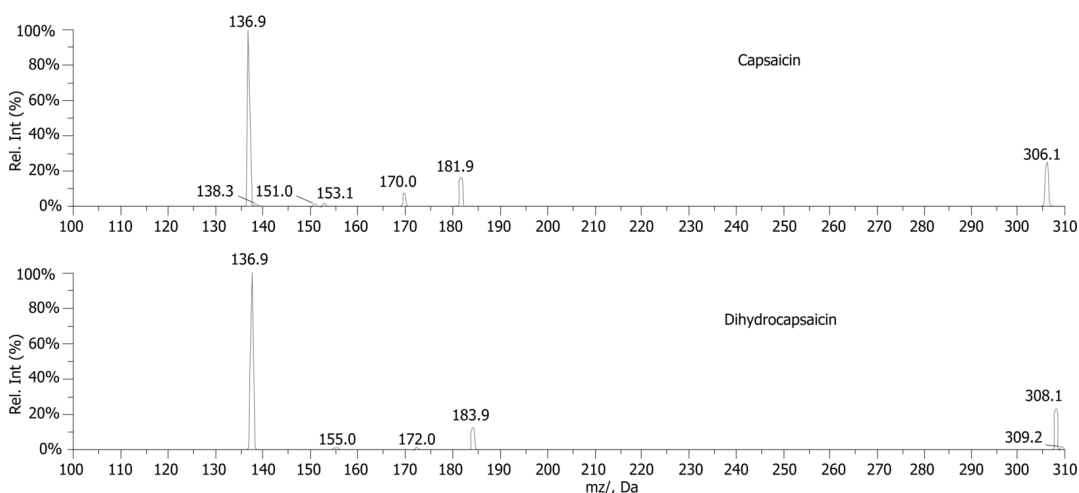

ANALYTES:

- 1 Capsaicin
- 2 Dihydrocapsaicin



Products used in this application:

