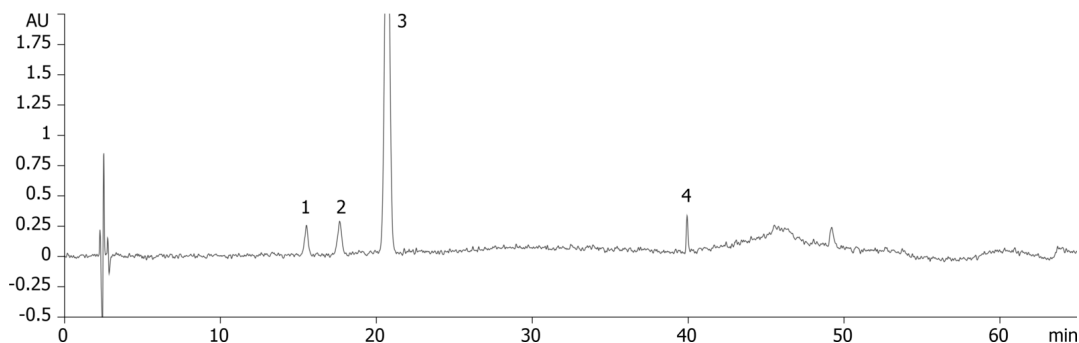


Paroxetine EP method on Kinetex 5 μ m C8 250x4.6 mm**Column:** Kinetex® 5 μ m C8 100 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4608-E0**Elution Type:** Gradient**Eluent A:** TFA/THF/Water (5:100:900)**Eluent B:** TFA/THF/ACN (5:100:900)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	30	80	20
	3	50	20	80
	4	55	20	80
	5	60	80	20
	6	65	80	20

Flow Rate: 1 mL/min**Col. Temp.:** 40 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

21927

**Products used in this application:****ANALYTES:**

- 1 Impurity A
- 2 Impurity E
- 3 Paroxetine
- 4 Impurity C

