

Pain Panel Opiates using Kinetex 2.6u Biphenyl 50x3.0mm

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4622-Y0

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Methanol with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	0.5	90	10
	3	2	75	25
	4	4.5	20	80
	5	4.51	15	85
	6	5.5	15	85
	7	7	90	10

Flow Rate: 0.6 mL/min


Col. Temp.: 40 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (40 °C)

22108



Products used in this application:

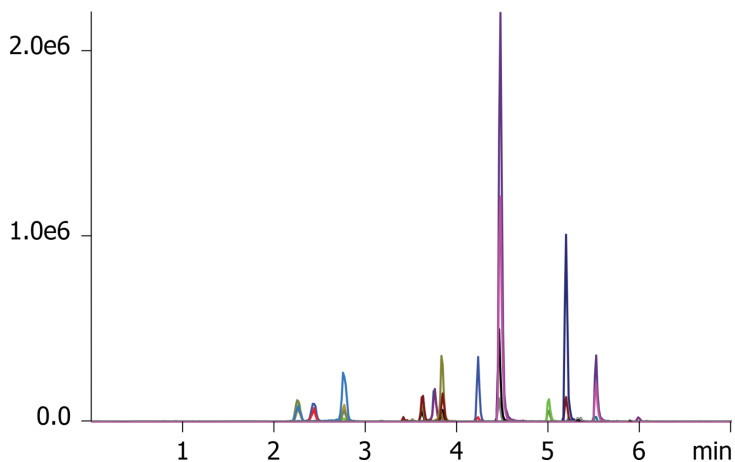

Pain Panel Opiates using Kinetex 2.6u Biphenyl 50x3.0mm

ANALYTES:

- 1** Morphine
- 2** Oxymorphone
- 3** Hydromorphone
- 4** Naloxone
- 5** Codeine
- 6** 6-MAM
- 7** Oxycodone
- 8** Naltrexone
- 9** Hydrocodone
- 10** Norfentanyl
- 11** Tramadol
- 12** Normeperidine
- 13** Meperidine
- 14** Norbuprenorphine
- 15** Fentanyl
- 16** Buprenorphine
- 17** Propoxyphene
- 18** Norpropoxyphene
- 19** EDDP
- 20** Methadone

