

Yarra 3µm SEC-2000 Flow Rate Study van Deemter curve


Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4512-K0
Elution Type: Isocratic
Eluent A: 100 mM Sodium Phosphate pH 6.8 + 0.02% NaN₃
Gradient
Step No. **Time (min)** **Pct A**
Profile: **1** 0 100
Flow Rate: 0.2 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
22126



Products used in this application:



Flow Rate vs. Proteins Efficiency on Yarra 2000

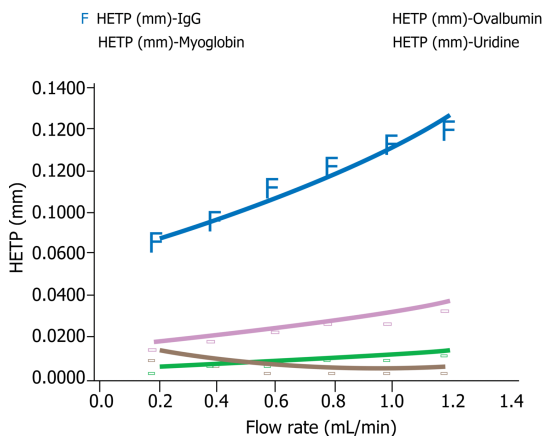

ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine

