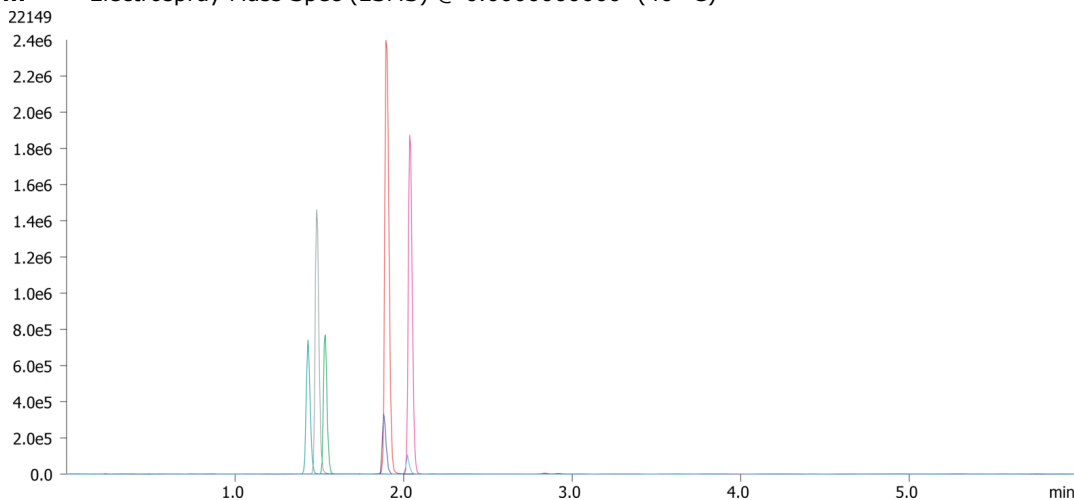


Antiepileptic Drugs by LCMS using Kinetex 2.6u Biphenyl 50x2.1mm**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** Water with 0.1% formic acid**Eluent B:** Acetonitrile with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	3	5	95
	3	3.1	98	2
	4	6	98	2

Flow Rate: 0.5 mL/min**Col. Temp.:** 40 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (40 °C)**ANALYTES:**

- 1** Pregabalin
Retention Time: 1.43 min
- 2** Levetiracetam
Retention Time: 1.48 min
- 3** Gabapentin
Retention Time: 1.53 min
- 4** Lamotrigine
Retention Time: 1.88 min
- 5** Lacosamide
Retention Time: 1.9 min
- 6** Zonisamide
Retention Time: 2.02 min
- 7** MHD
Retention Time: 2.04 min

**Products used in this application:**