

TCA's using Kinetex 5u C18 150x4.6mm


Column: Kinetex® 5µm C18 100 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4601-E0
Elution Type: Isocratic
Eluent A: 20mM Potassium Phosphate pH7.0
Eluent B: 50:50 Acetonitrile:Methanol
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	30	70

Profile:
Flow Rate: 1.5 mL/min
Col. Temp.: 30 °C
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)



Products used in this application:

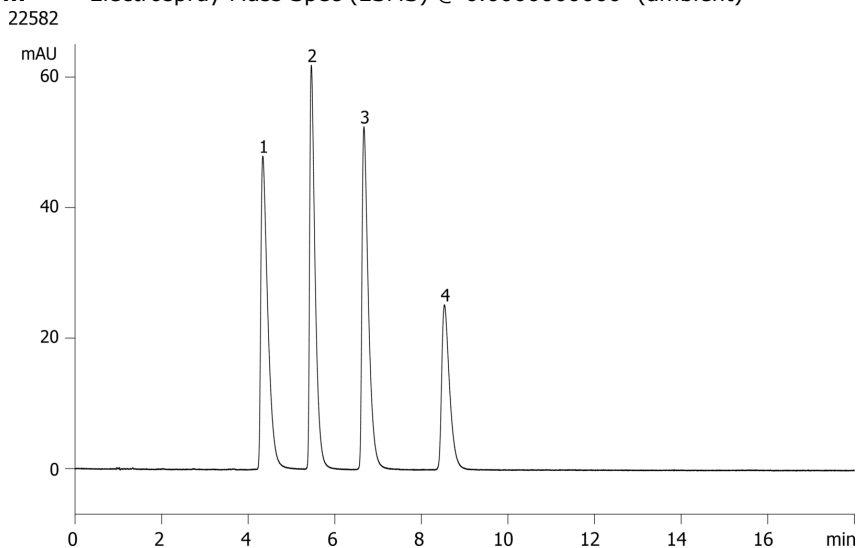

ANALYTES:

- 1** Imipramine
Retention Time: 4.337 min
- 2** Nortriptyline
Retention Time: 5.458 min
- 3** Amitriptyline
Retention Time: 6.673 min
- 4** Clomipramine
Retention Time: 8.533 min

