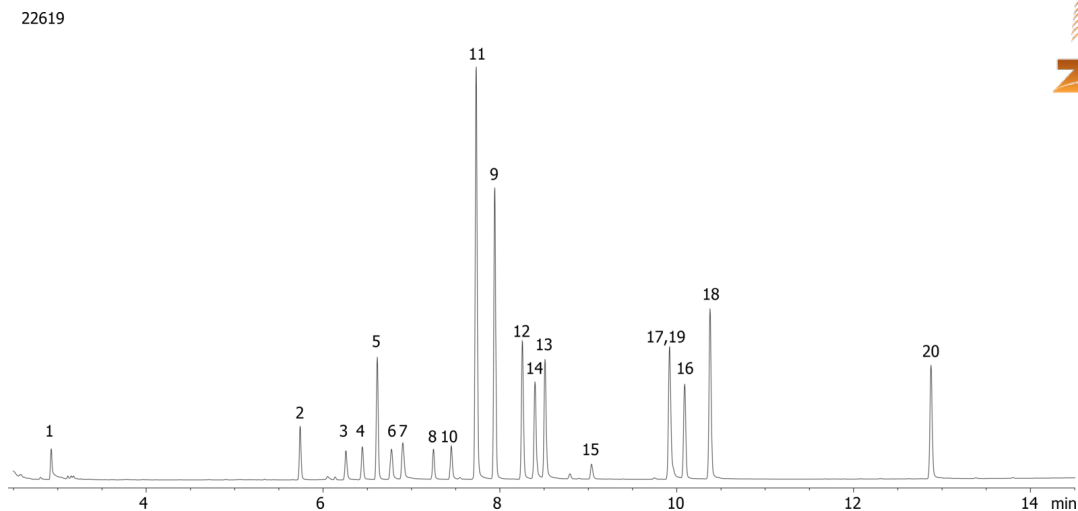


EPA 8151 by GC/ECD on ZB-CLPesticides-1 (Helium Gas)

Column: Zebtron ZB-CLPesticides-1 Column 30m x 0.32mm x 0.32 μ m , Ea
Phase: Proprietary Pesticides Phase
Dimensions: 30 meters x 0.32 mm x 0.32 μ m
Order No: 7HM-G028-51
Oven Profile: 70 $^{\circ}$ C for 0.5 min to 190 $^{\circ}$ C @ 25 $^{\circ}$ C/min for 1 min to 300 $^{\circ}$ C @ 11 $^{\circ}$ C/min for 5 min
Carrier Gas: Constant Flow Helium, 36 cm/sec
Injection: Splitless 0.000000000:1 0 μ L @ 250 $^{\circ}$ C
Detection: Refractive Index (325 $^{\circ}$ C)
Analyst Note: Analytes are 0.1 ug/mL in Hexane

**Products used in this application:**

EPA 8151 by GC/ECD on ZB-CLPesticides-1 (Helium Gas)

ANALYTES:

- 1 Dalapon methyl ester
- 2 3,5-Dichlorobenzoic acid methyl ester (SS)
- 3 4-Nitroanisole
- 4 DCAA methyl ester (SS)
- 5 Dicamba methyl ester
- 6 MCPP methyl ester
- 7 MCPA Methyl ester
- 8 Dichlorprop, methyl ester
- 9 4,4'-DBOB (IS)
- 10 2,4-D methyl ester
- 11 Pentachloroanisole
- 12 2,4,5-TP, methyl ester
- 13 2,4,5-T methyl ester
- 14 Chloramben, methyl ester
- 15 2,4-DB methyl ester
- 16 Dinoseb methyl ester
- 17 Bentazon methyl ester
- 18 DCPA methyl ester (Chlorthal-dimethyl)
- 19 Pichloram methyl ester
- 20 Acifluorfen methyl ester

