

Pharma ABN Mixture using Kinetex 5u EVO C18 150x4.6mm

Column: Kinetex® 5µm EVO C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4633-E0

Elution Type: Gradient

Eluent A: Water with 0.1% Formic acid


Eluent B: Acetonitrile with 0.1% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	5	95
	3	10.1	95	5
	4	15	95	5

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)

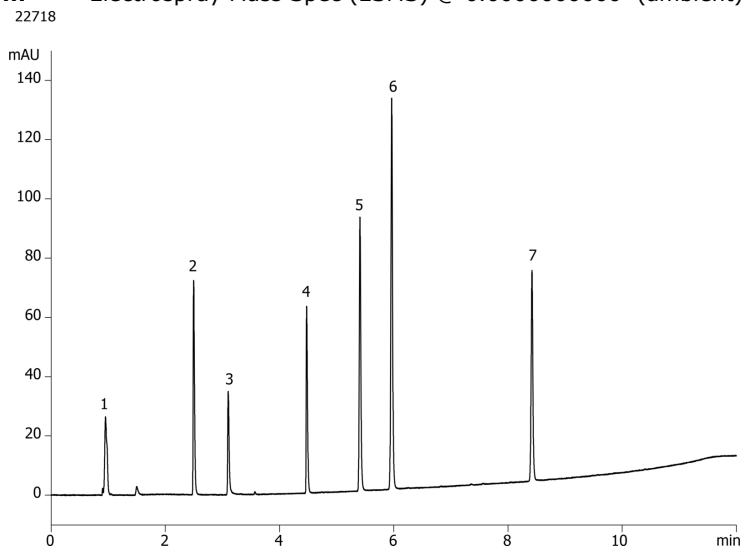

ANALYTES:

- 1** Uracil
Retention Time: 0.953 min
- 2** Pindolol
Retention Time: 2.497 min
- 3** Chlorpheniramine
Retention Time: 3.102 min
- 4** Nortriptyline
Retention Time: 4.475 min
- 5** 3-Methyl-4-Nitrobenzoic acid
Retention Time: 5.409 min
- 6** 5-Methyl Salicyl Aldehyde
Retention Time: 5.964 min
- 7** Hexanophenone
Retention Time: 8.421 min



Products used in this application:

