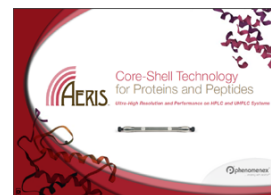
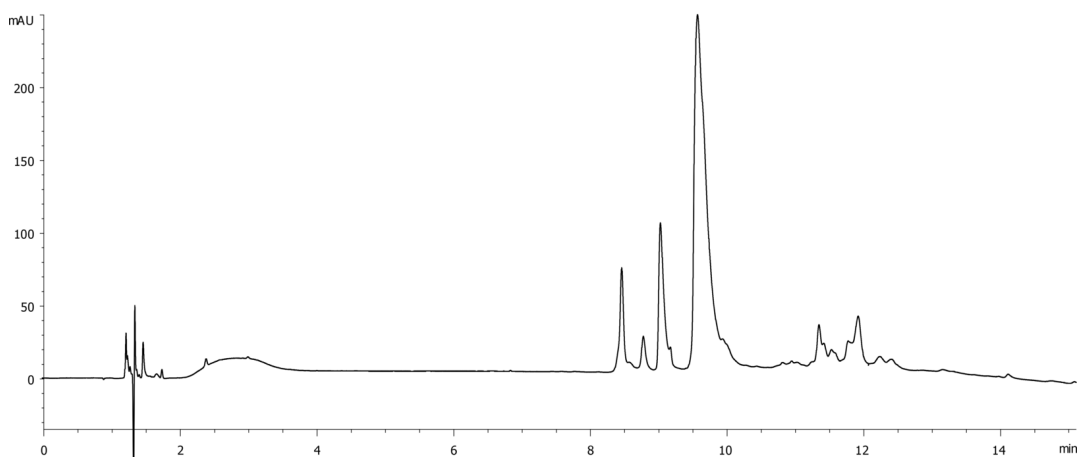


Aeris Ig-G comparison widepore XB-C8**Column:** Aeris[™] 3.6 μ m WIDEPORE XB-C8 200 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4481-E0**Elution Type:** Gradient**Eluent A:** 0.1% TFA**Eluent B:** 0.085% TFA in AcCN

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	2	1	80	20
	3	15	40	60
	4	15.01	90	10

Flow Rate: 1 mL/min**Col. Temp.:** 40 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

22746

**Products used in this application:****ANALYTES:****1** Ig-G1 Standard