

## Antibody Fluorophore Conjugate (FITC-IgG) Low Salt Buffer on Yarra 3u SEC-3000 300 x 7.8 mm


**Column:** Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea  
**Dimensions:** 300 x 7.8 mm ID  
**Order No:** 00H-4513-K0  
**Elution Type:** Isocratic  
**Eluent A:** 100mM PO4 pH 6.8  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)  
22839



Products used in this application:

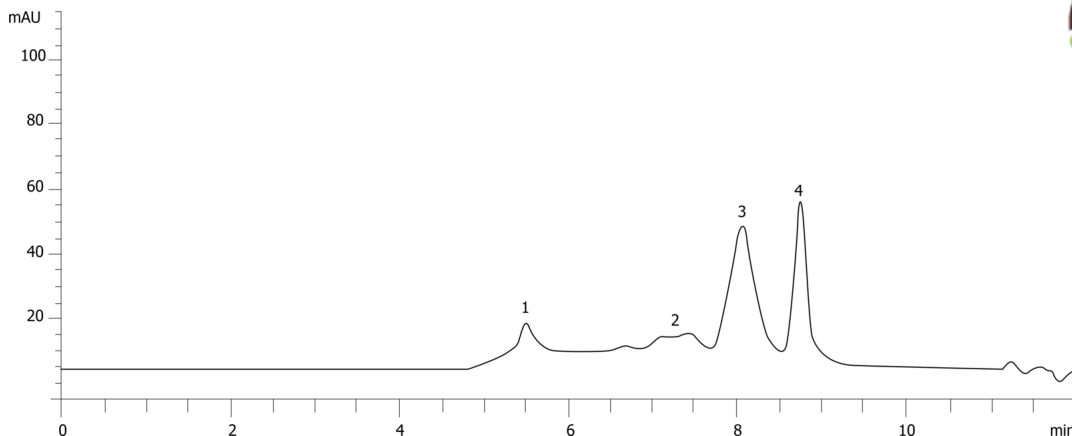

### ANALYTES:

- 1 HMW Aggregate
- 2 Dimer
- 3 FITC-IgG
- 4 LMW Fragment

