

Caffeine in Dietary Supplement Tablets by HPLC

Column: Synergi™ 4 µm Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4337-E0

Elution Type: Isocratic

Eluent A: Acetonitrile/50 mM Sodium Acetate, 10:90

Gradient


Profile:

Step No.	Time (min)	Pct A
1	10	100

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (25 °C)

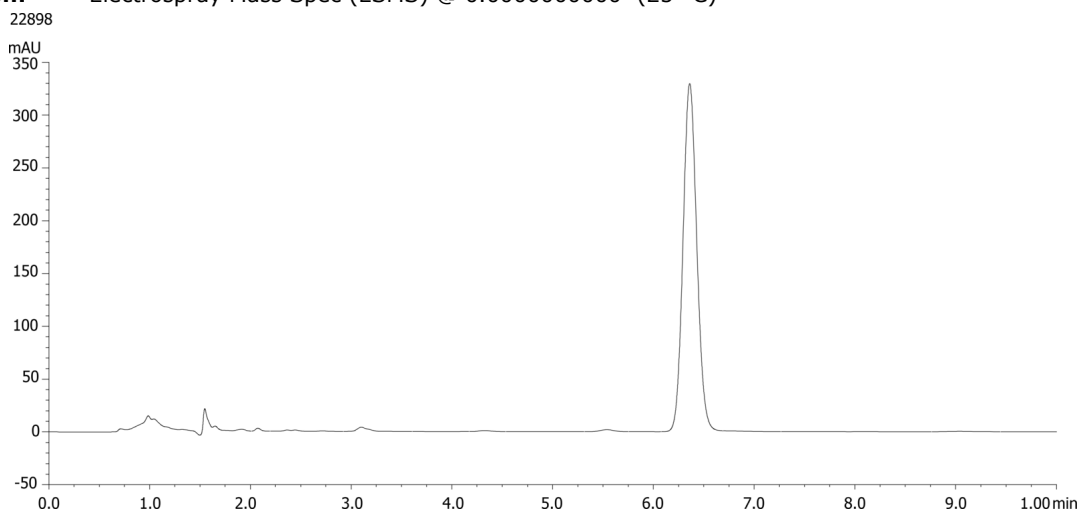

ANALYTES:

1 Caffeine

Retention Time: 6.5 min

