

Vitamin E Acetate in Dietary Supplement Tablets by HPLC

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile/Methanol, 75:25


Gradient	Step No.	Time (min)	Pct A
----------	----------	------------	-------

Profile:	1	20	100
-----------------	----------	----	-----

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (25 °C)

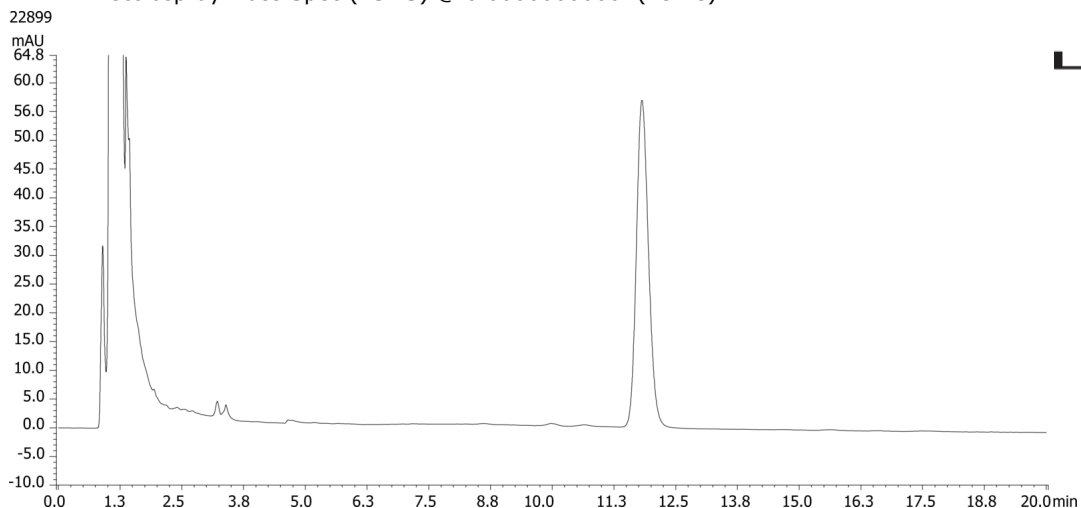

Products used in this application:

**ANALYTES:**

1 Vitamine E Acetate

Retention Time: 13.6 min