

Cold medicine compounds on Kinetex 2.6 µm EVO C18

Column: Kinetex® 2.6µm EVO C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4725-E0

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.1% TFA in ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	6	70	30
	3	11	0	100

Flow Rate: 1.25 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

Analyst Note: Forced degradation protocol:


Reconstitute 1 mg of cefaclor neat standard with 1.0 mL of 1 M HCl. Heat the solution at 60 °C for 6 hours then neutralize with 1.0 mL of 1 M NaOH. Refrigerate until ready to inject.



Products used in this application:



23036

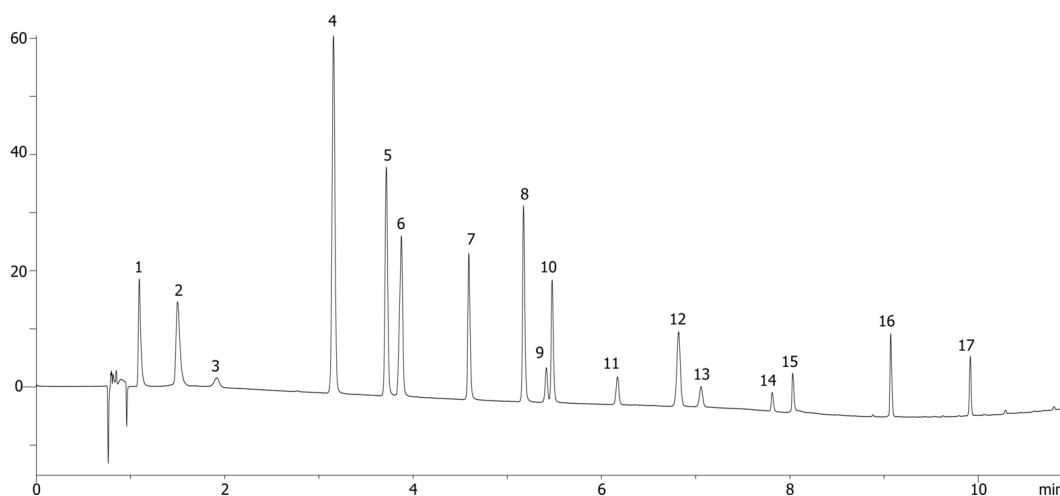

Cold medicine compounds on Kinetex 2.6 µm EVO C18

ANALYTES:

- 1** Maleic acid
- 2** Fumaric acid
- 3** Phenylephrine
- 4** Acetaminophen
- 5** Pheniramine
- 6** Doxylamine
- 7** Pyrilamine
- 8** Chlorpheniramine
- 9** Guaifenesin
- 10** Brompheniramine
- 11** Acetylsalicylic acid
- 12** 4-Nitrophenol
- 13** Impurity
- 14** Dextromethorphan
- 15** Diphenhydramine
- 16** Clemastine
- 17** Ibuprofen

