

Fast Analysis of Sugars in Fruit Juices - Moscato (White Wine)

Column: Rezex™ RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

Dimensions: 100 x 7.8 mm ID

Order No: 00D-0130-K0

Elution Type: Isocratic

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.4 mL/min

Col. Temp.: 80 °C


Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (35 °C)
23056

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

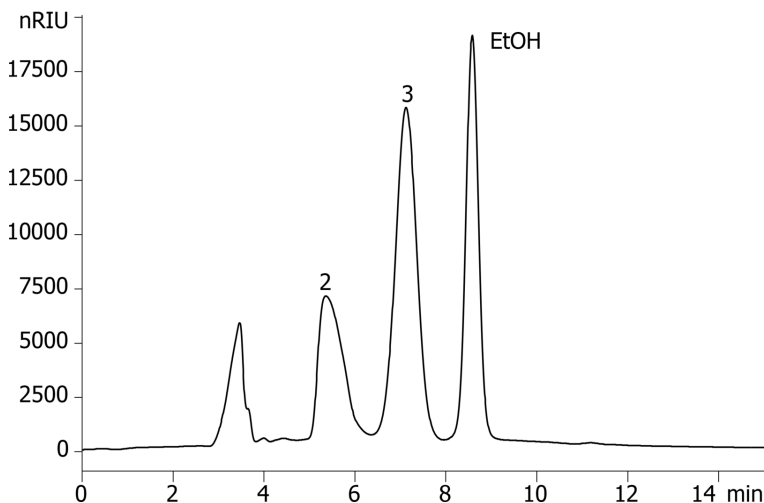

ANALYTES:

1 Moscato (White Wine)

