


Pesticides in Poultry Meat by GC/ECD on ZB-MultiResidue-1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G016-17**Oven Profile:** 150 $^{\circ}$ C for 0.5 min to 20 $^{\circ}$ C/min to 220 $^{\circ}$ C hold 0 min, 6 $^{\circ}$ C, to 310 $^{\circ}$ C hold time: 0 min**Carrier Gas:** Ramp Flow Helium, 3.4 mL/min**Injection:** Splitless 0.0000000000:1 0 μ L @ 250 $^{\circ}$ C**Detection:** Refractive Index (350 $^{\circ}$ C)**Analyst Note:** Flow ramp program: 3.4mL/min ramped flow- Hold 6.00 min, ramp 5.00 mL/min to 8.00 mL/min, hold 12.08 min,**Products used in this application:**

Pesticides in Poultry Meat by GC/ECD on ZB-MultiResidue-1

ANALYTES:

- 1** 2,4,5,6-tetrachloro-m-xylene (Internal Standard 1)
- 2** aBHC
- 3** HCB
- 4** Lindane
- 5** Heptachlor
- 6** Drusban
- 7** Aldrin
- 8** Heptachlor Epoxide
- 9** p,p'-DDE
- 10** Dieldrin
- 11** o,p-DDT
- 12** p,p'-TDE
- 13** Endrin
- 14** p,p'-DDT
- 15** Methoxychlor
- 16** Merex
- 17** Decachlorobiphenyl (Internal Standard 2)

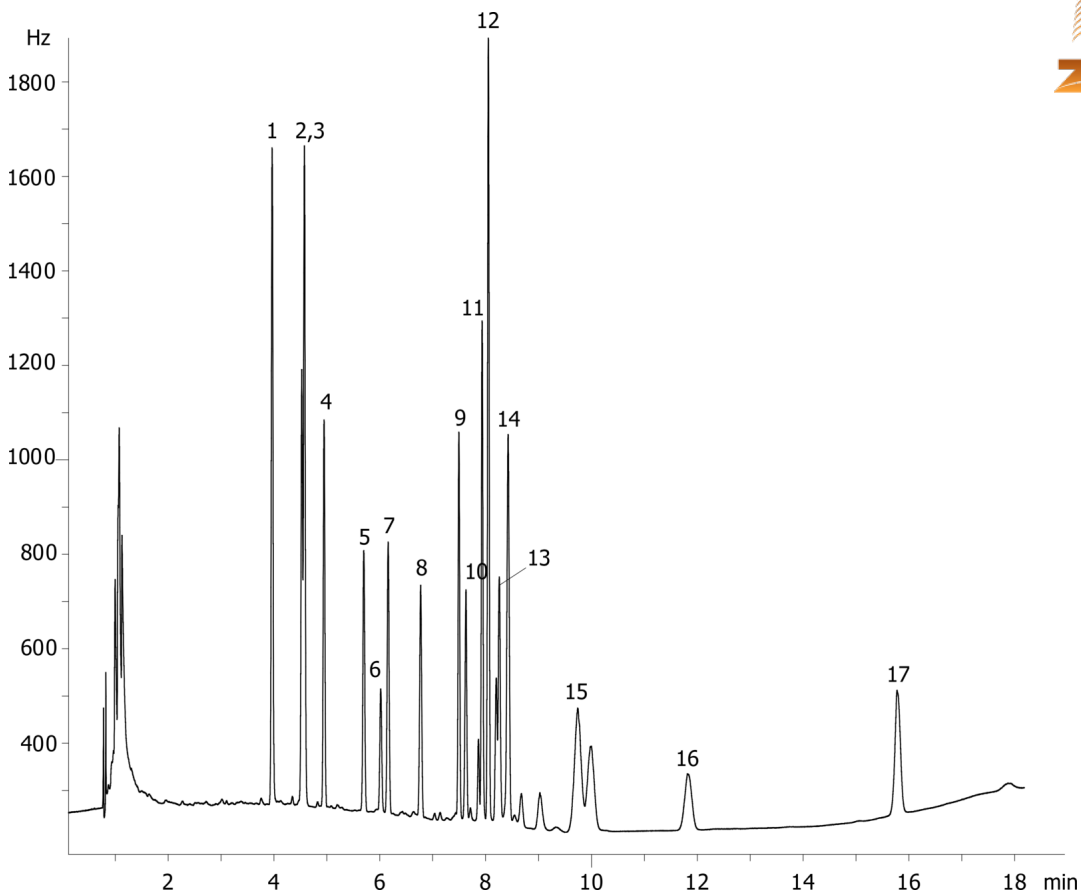