

Mirtazapine on LUX Amylose-1 3u 250x4.6mm in NP

Column: Lux® 3 µm Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4729-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% DEA in IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	10	80	20


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)



Products used in this application:

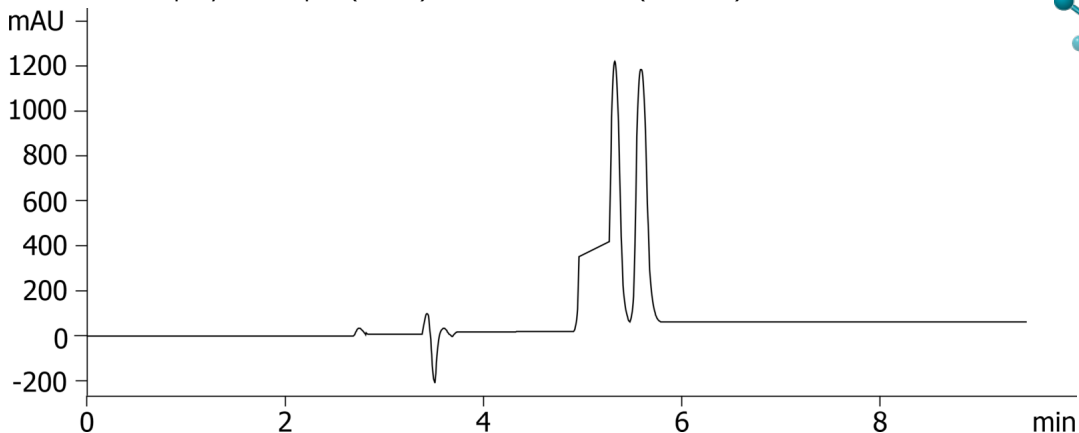

ANALYTES:

- 1 Mirtazapine
Retention Time: 5.319 min
- 2 Mirtazapine
Retention Time: 5.564 min