

Drugs of Abuse by LCMS using Kinetex 1.7u Biphenyl 50x2.1mm

Column: Kinetex® 1.7 µm Biphenyl 100 Å LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4628-AN

Elution Type: Gradient

Eluent A: 0.1 Formic Acid in Water


Eluent B: 0.1 Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	5	95
	3	5	5	95
	4	5.1	95	5

Flow Rate: 0.4 mL/min

Col. Temp.: 40 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

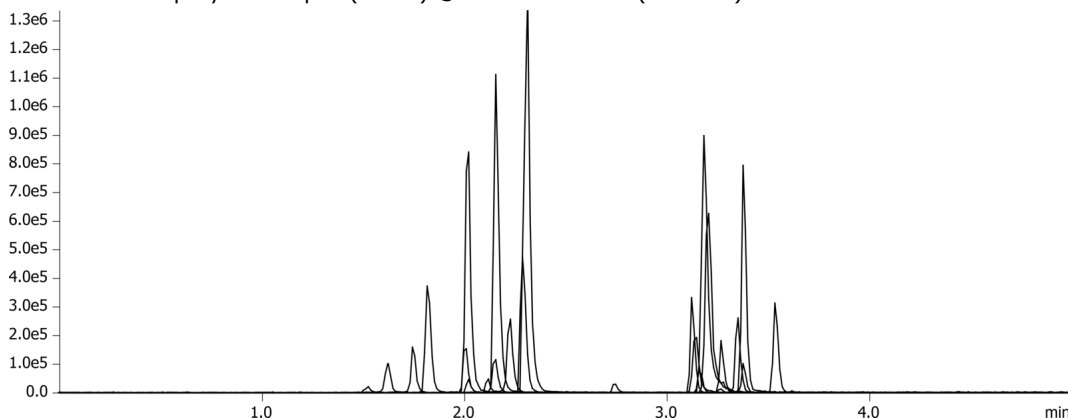

Products used in this application:



Drugs of Abuse by LCMS using Kinetex 1.7u Biphenyl 50x2.1mm

ANALYTES:

- 1** Morphine
- 2** Oxymorphone
- 3** Hydromorphone
- 4** Amphetamine
- 5** MDA
- 6** Codeine
- 7** Methamphetamine
- 8** 6-MAM
- 9** Oxycodone
- 10** MDMA
- 11** Hydrocodone
- 12** BZE
- 13** MDEA
- 14** Norbuprenorphine
- 15** Oxazepam
- 16** Midazolam
- 17** Lorazepam
- 18** PCP
- 19** Flurazepam
- 20** Buprenorphine
- 21** Clonazepam
- 22** Alprazolam
- 23** Temazepam
- 24** Diazepam
- 25** Methadone

