


Pyrethroids Extracted from Marine Sediment using QuEChERS by LCMSMS**Column:** Kinetex® 2.6µm EVO C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4725-AN**Elution Type:** Gradient**Eluent A:** 5mM ammonium acetate**Eluent B:** methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50
	2	1	5	95
	3	3	5	95
	4	3.01	50	50
	5	6	50	50

Flow Rate: 700**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)**Products used in this application:****ANALYTES:**

- 1 Cyfluthrin
- 2 Cyhalothrin
- 3 Cypermethrin
- 4 Deltamethrin
- 5 Esfenvalerate
- 6 cis-Permethrin
- 7 trans-Permethrin
- 8 Bifenthrin

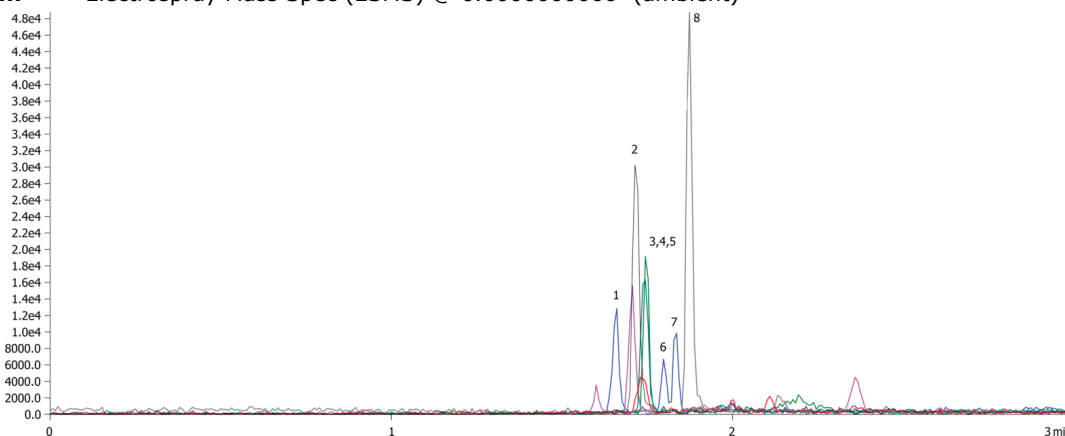