

Amitriptyline gradient loading using Luna Omega 1.6u PS C18 50x2.1mm

Column: Luna Omega 1.6u PS C18 100A, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4752-AN

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Acetonitrile with 0.1% formic acid

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------


Profile: 1 0 95 5

2 5 20 80


Flow Rate: 0.4 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

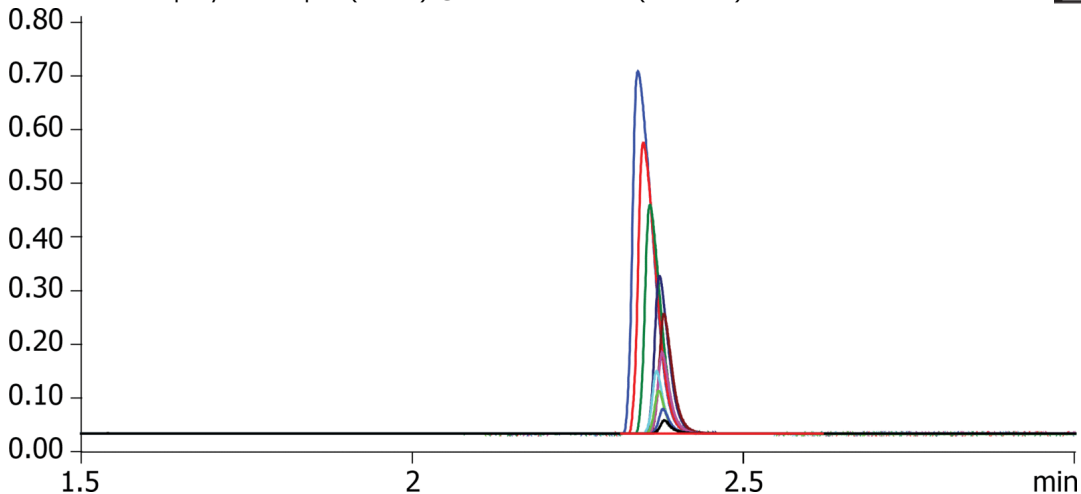**ANALYTES:**

1 Amitriptyline

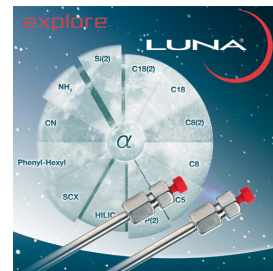

Products used in this application:

