

Chiral Separation of Chlorthalidone on Lux 5u i-Cellulose-5


Column: Lux® 5 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4756-E0
Elution Type: Isocratic
Eluent A: 0.1% Ethylamine in MtBE
Eluent B: 0.1% Ethylamine in MeOH
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	96	4

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)
 23758



Products used in this application:

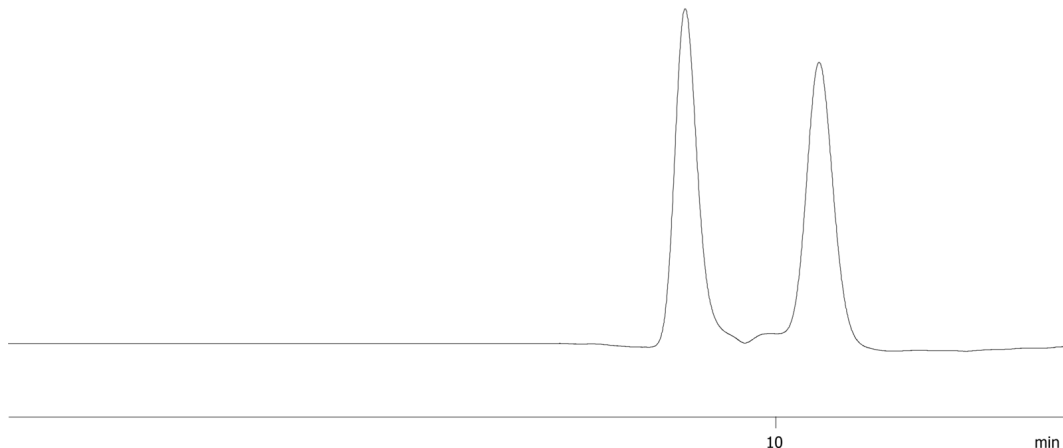

ANALYTES:

- 1 Chlorthalidone
Retention Time: 8.787 min
- 2 Chlorthalidone
Retention Time: 10.564 min

