

Chiral Separation of Althiazide on Lux 3u i-Cellulose-5 (Batch#A)


Column: Lux® 3 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4755-E0
Elution Type: Isocratic
Eluent A: 0.1% Ethylamine in Dichloromethane
Eluent B: 0.1% Ethylamine in MeOH
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	98	2

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
23764



Products used in this application:

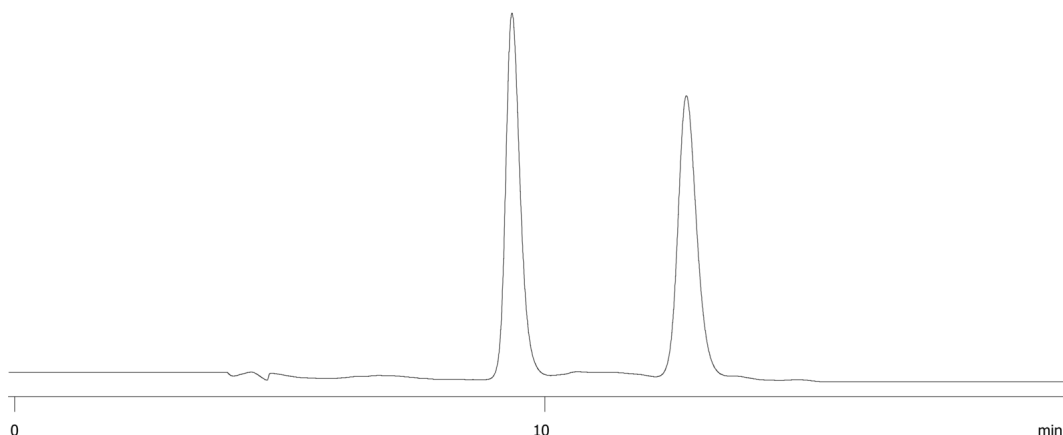

ANALYTES:

- 1 Althiazide
Retention Time: 9.483 min
- 2 Althiazide
Retention Time: 12.87 min

