

USP system suitability for Atomoxetine HCl and Amoxetine S Isomer on Lux i-Cellulose-5


Column: Lux® 5 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4756-E0
Elution Type: Isocratic
Eluent A: 970:30:5 Hexane:IPA:DEA
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
23768



Products used in this application:

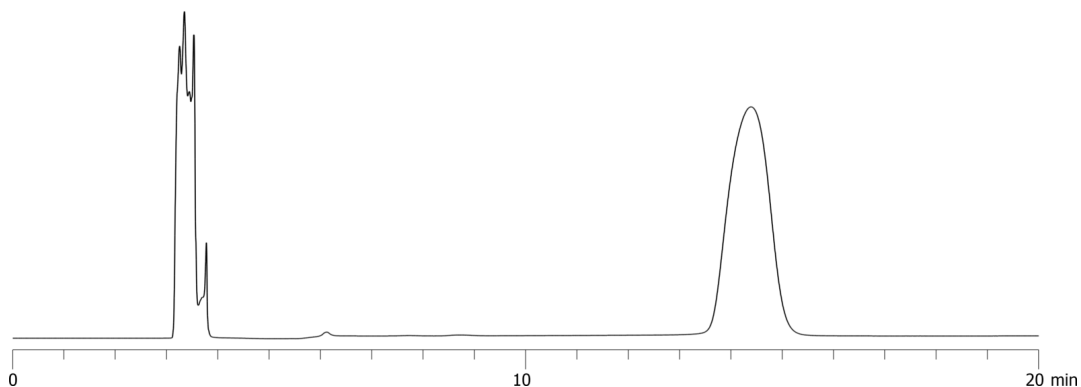

ANALYTES:

- 1 Atomoxetine S Isomer
Retention Time: 8.1 min
- 2 Atomoxetine HCl
Retention Time: 11.4 min

