

FAMES In Infant Formula by GC/FID on Zebron ZB-FAME

Column: Zebron ZB-FAME, GC Cap. Column 30m x 0.25mm x 0.2µm , Ea
Phase: Proprietary Pesticides Phase
Dimensions: 30 meters x 0.25 mm x 0.2 µm
Order No: 7HG-G033-10
Oven Profile: 100 Å°C for 2 min to 140 Å°C @ 10 Å°C/min to 190 Å°C @ 3 Å°C/min to 260 Å°C @ 30 Å°C/min for 2 min

Carrier Gas: Constant Flow Helium, 1.2 mL/min
Injection: Split 50.0000000000:1 0 µL @ 240°C
Detection: Refractive Index (260°C)

Analyst Note: Recommended Liner: Zebron PLUS Single Taper with Wool, 4 mm ID
Liner Part No.: AG2-0A11-05 (for Agilent systems)
Inlet Seal: AG0-8620 (Gold Plated Easy Seal)
Septum: AG0-4696 (PhenoRed-400)

Sample Preparation:


1. Weigh out 500 mg of powdered infant formula into a scintillation vial
2. Dissolve (or dilute) sample in 5 mL toluene
3. Add 6 mL 10 % acetyl chloride solution in methanol
4. Incubate @ 80 °C for 2 hours
5. After incubation add 10 mL 10 % Na₂CO₃ solution and centrifuge at 5000 rpm for 5 min
6. Extract organic layer for GC analysis



Products used in this application:



23827

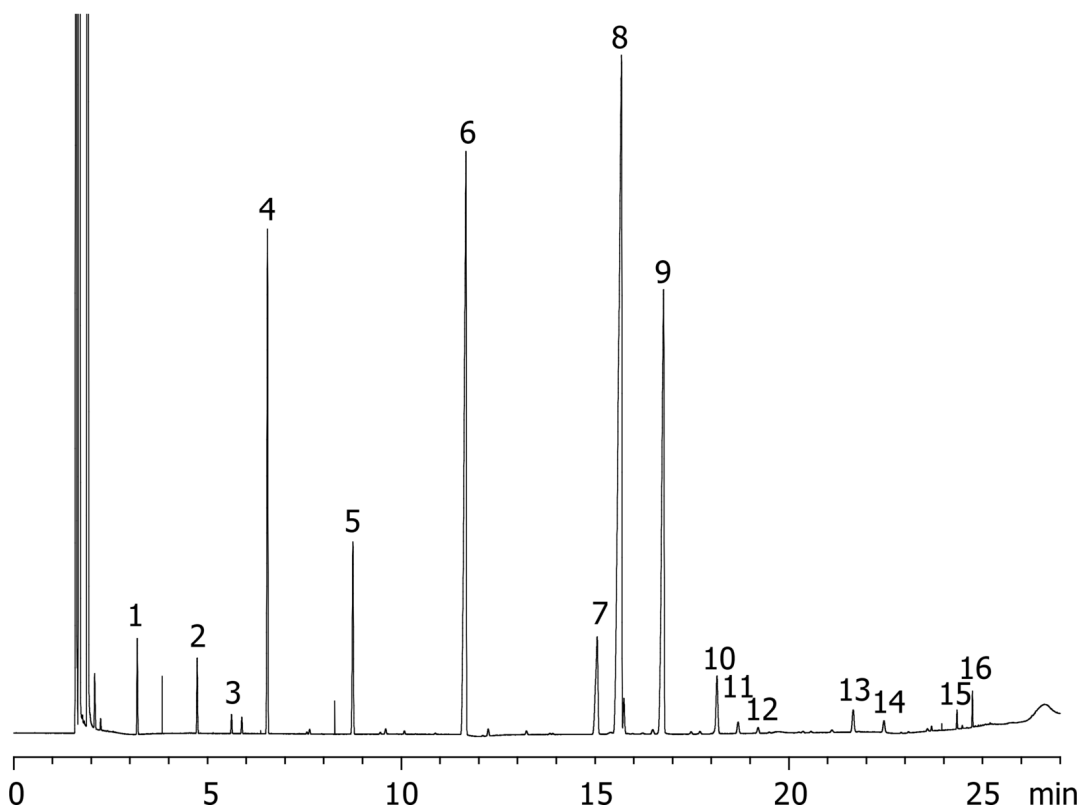

FAMES In Infant Formula by GC/FID on Zebron ZB-FAME

ANALYTES:

- 1 C8:0
- 2 C10:0
- 3 C11:0 (I.S.)
- 4 C12:0
- 5 C14:0
- 6 C16:0
- 7 C18:0
- 8 C18:1 cis 9
- 9 C18:2 cis 9,12
- 10 C18:3 cis 9,12,15
- 11 C20:0
- 12 C20:1 cis 11
- 13 C20:4 cis 5,8,11,14
- 14 C22:0
- 15 C24:0
- 16 C22:6 cis 4,7,10,13,16,19

