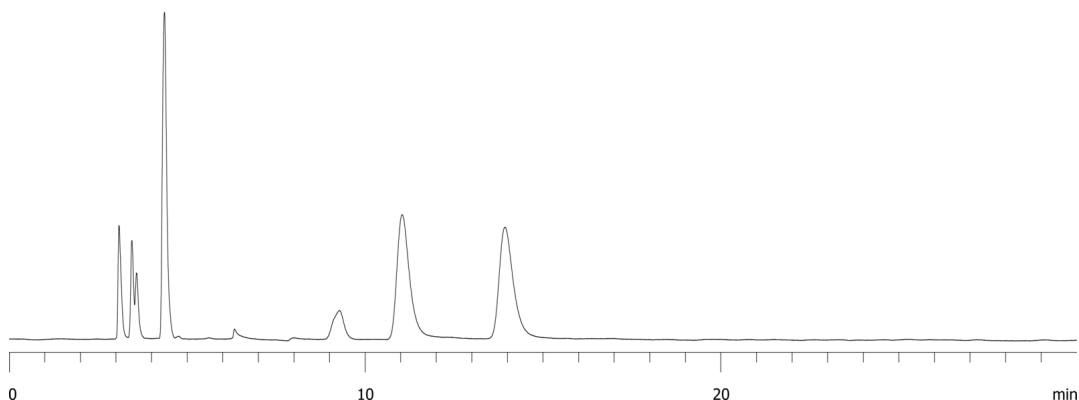


Chiral Separation of 2-Chloromandelic Acid on Lux 5 μ m i-Cellulose-5**Column:** Lux[®] 5 μ m i-Cellulose-5, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4756-E0**Elution Type:** Isocratic**Eluent A:** 0.1% TFA in Dichloromethane**Eluent B:** 0.1% TFA in IPA**Gradient**

Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	99	1

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)

23879

**Products used in this application:****ANALYTES:****1** 2-Chloromandelic acid

Retention Time: 11 min

2 2-Chloromandelic acid

Retention Time: 13.9 min

