

Chiral Separation of Tropic Acid on Lux 3µm i-Cellulose-5


Column: Lux® 3 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4755-E0
Elution Type: Isocratic
Eluent A: 0.1% TFA in Hexane
Eluent B: 0.1% TFA in Ethyl acetate
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	70	30

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
 23882



Products used in this application:

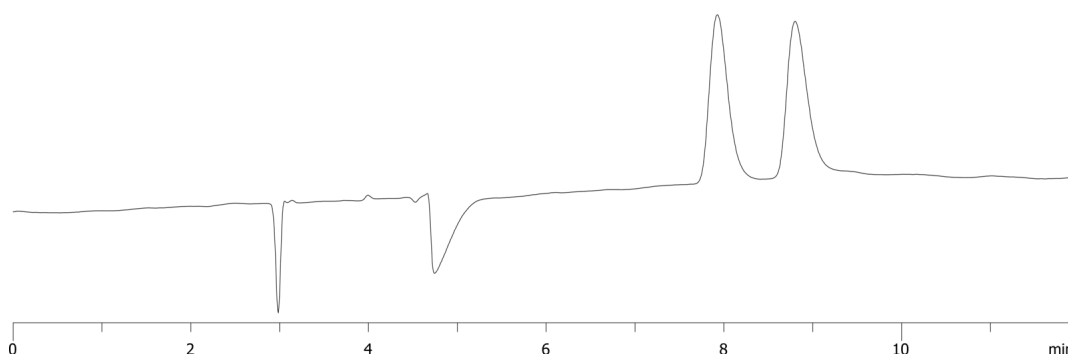

ANALYTES:

- 1 Tropic acid
Retention Time: 7.93 min
- 2 Tropic acid
Retention Time: 8.8 min

