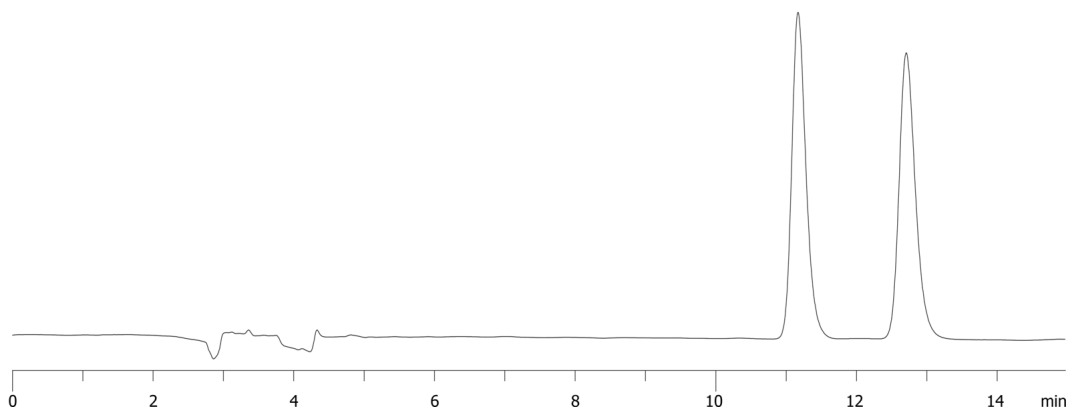


Chiral Separation of Tetramisole on Lux 3 μ m i-Cellulose-5

Column: Lux[®] 3 μ m i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4755-E0
Elution Type: Isocratic
Eluent A: 0.1% DEA in Hexane
Eluent B: 0.1% DEA in Ethanol
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	70	30

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
23886

**Products used in this application:****ANALYTES:**

- 1 Tetramisole
Retention Time: 11.16 min
- 2 Tetramisole
Retention Time: 12.7 min

