

Chiral Separation of Sulpiride on Lux 3µm i-Cellulose-5

Column: Lux® 3 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4755-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% DEA in THF

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	60	40

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)

23888



Products used in this application:

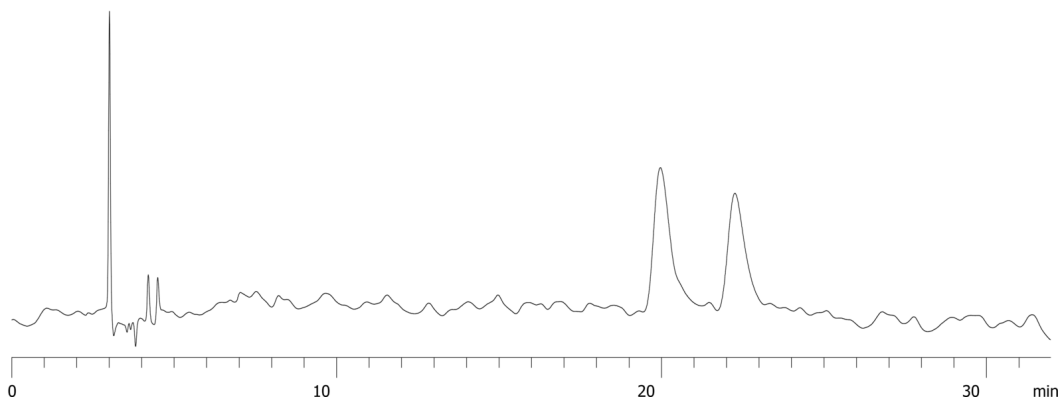

ANALYTES:

1 Sulpiride

Retention Time: 19.96 min

2 Sulpiride

Retention Time: 22.254 min

