

Chiral Separation of Nicardipine on Lux 5µm i-Cellulose-5


Column: Lux® 5 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4756-E0
Elution Type: Isocratic
Eluent A: 0.1% DEA in Hexane
Eluent B: 0.1% DEA in IPA
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	90	10

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)
 23891



Products used in this application:

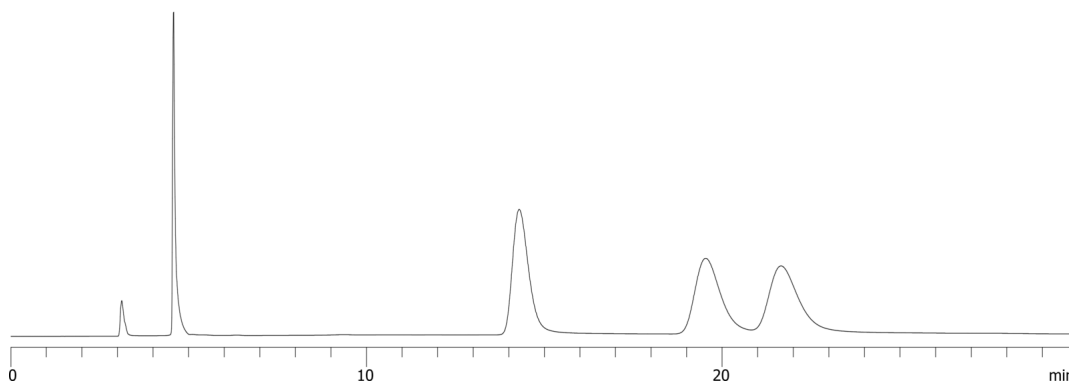

ANALYTES:

- 1** Nicardipine
Retention Time: 14.289 min
- 2** Nicardipine
Retention Time: 19.523 min
- 3** Nicardipine
Retention Time: 21.65 min

