

FAMES In Canola Oil by GC/MS on Zebron ZB-FAME

Column: Zebron ZB-FAME, GC Cap. Column 30m x 0.25mm x 0.2um , Ea

Phase:

Dimensions: 30 meters x 0.25 mm x 0.2 µm

Order No: 7HG-G033-10

Oven Profile: 100 Å°C for 2 min to 214 Å°C @ 3 Å°C/min

Carrier Gas: Constant Flow Helium, 1 mL/min

Injection: Split 50.0000000000:1 0 µL @ 240°C

Detection: Refractive Index (250°C)

Analyst Note: Recommended Accessories:
 Recommended Liner: Zebron PLUS Single Taper Z-Liner, 4 mm ID
 Liner Part No.: AG2-0A13-05 (for Agilent systems)
 Inlet Seal: AG0-8620 (Gold Plated Easy Seal)
 Septum: AG0-4696 (PhenoRed-400)

Sample Preparation:


1. Weigh out 25 mg of oil into a 5 mL screw cap glass vial
2. Add 2 mL of BF₃-methanol to sample vial
3. Heat vial at 60 °C for 15 minutes with periodic vortexing
4. Remove from heat and allow to attain room temperature
5. Add 1 mL water followed by 1 mL hexane
6. Vortex for 1 minute
7. Remove upper organic layer
8. Dry hexane layer over anhydrous sodium sulfate
9. Transfer to autosampler vial and inject on GC



Products used in this application:



24074

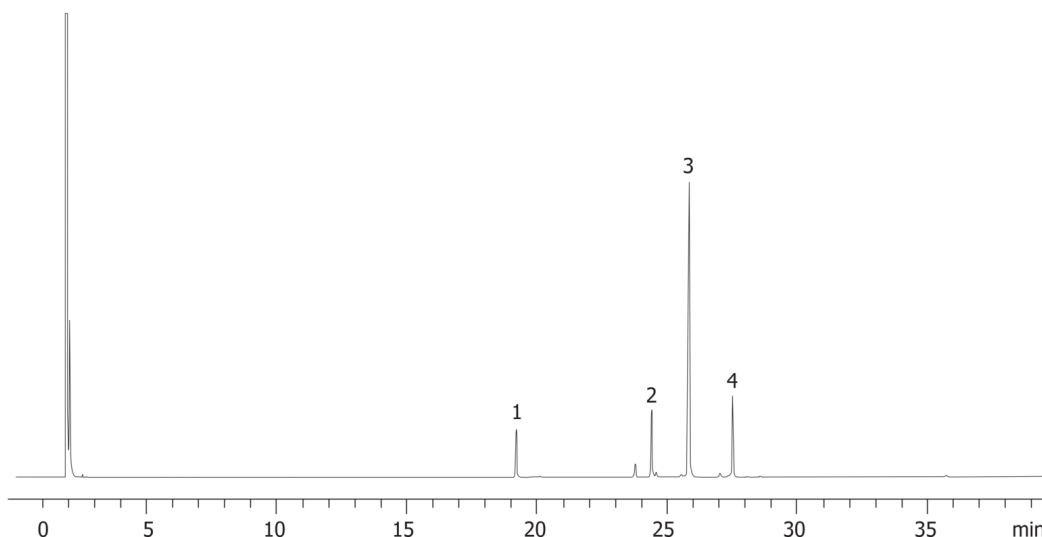

ANALYTES:

- 1 C16:0
- 2 C18:1 cis 9
- 3 C18:2 cis 9,12
- 4 C18:3 cis 9,12,15

©2015 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com

support@phenomenex.com