

Peptide analysis by LCMS using Luna Omega 3um PS C18 50x2.1mm

Column: Luna® Omega 3 µm PS C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4758-AN

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Acetonitrile with 0.1% formic acid

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	100	0
	2	4	0	100

Flow Rate: 0.5 mL/min


Col. Temp.: 25 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

24148



Products used in this application:

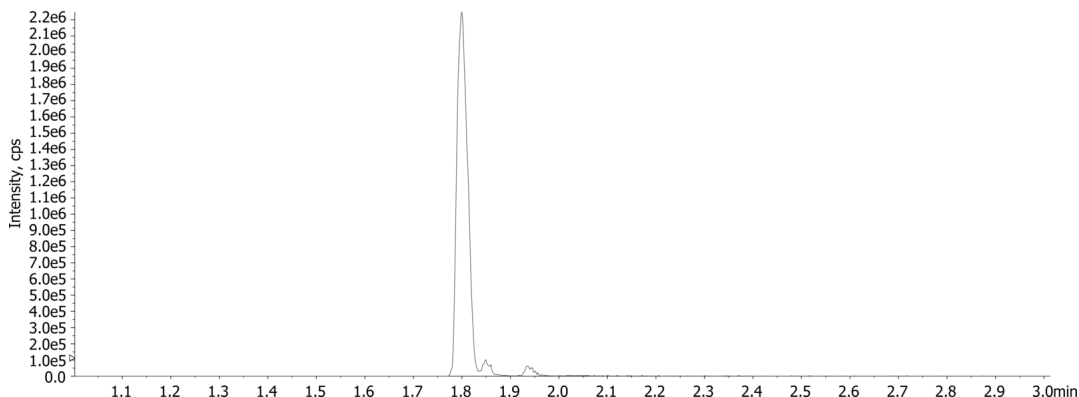

ANALYTES:

- 1 Angiotensin
- 2 Isobaric interference
- 3 Isobaric interference #2

