

Cromolyn sodium using Luna Omega 3um Polar C18 50x4.6mm

Column: Luna® Omega 3 µm Polar C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4760-E0

Elution Type: Gradient

Eluent A: 10mM ammonium formate pH 3.3

Eluent B: Acetonitrile with 0.1% formic acid

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	2	5	10	90

Flow Rate: 1.5 mL/min


Col. Temp.: 40 °C

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

24173



Products used in this application:

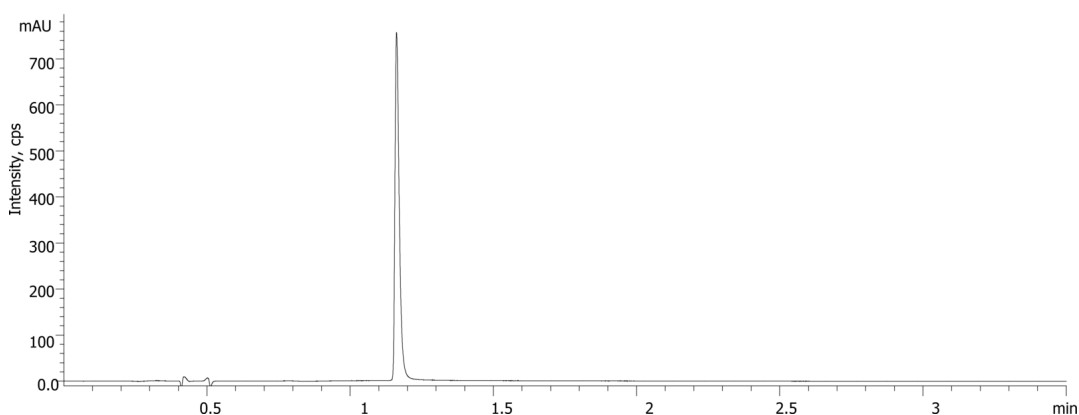

ANALYTES:

- 1 Cromolyn sodium (1mg/ml)

