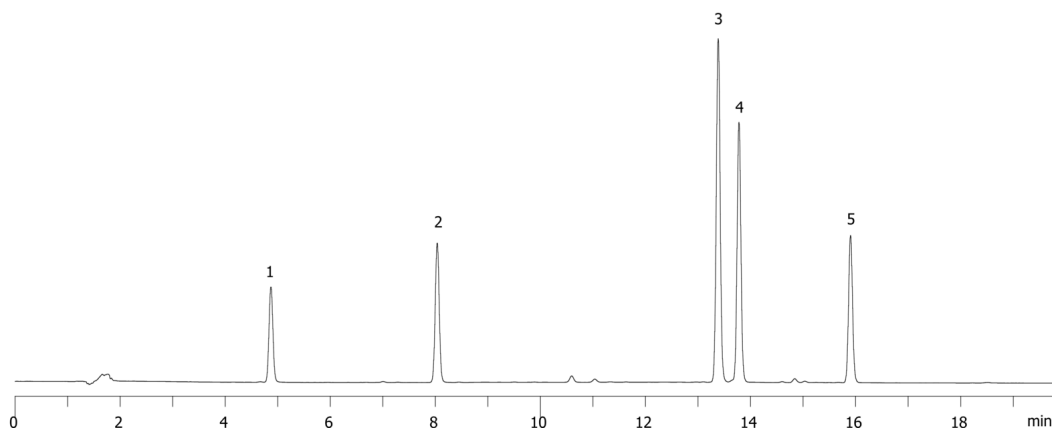


Taxanes on Kinetex 5 μ m F5**Column:** Kinetex® 5 μ m F5 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4724-E0**Elution Type:** Gradient**Eluent A:** Water + 0.1 % Formic Acid**Eluent B:** Acetonitrile + 0.1 % Formic Acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	20	30	70

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
24258**Products used in this application:****ANALYTES:****1** 7-Epi 10-Desacetyl Paclitaxel

Retention Time: 13.6 min

2 10-deacetylbaaccation III

Retention Time: 4.8 min

3 Baccatin III

Retention Time: 8 min

4 7-epi-taxol

Retention Time: 15.9 min

5 Taxol (Paclitaxel)

Retention Time: 13.9 min

