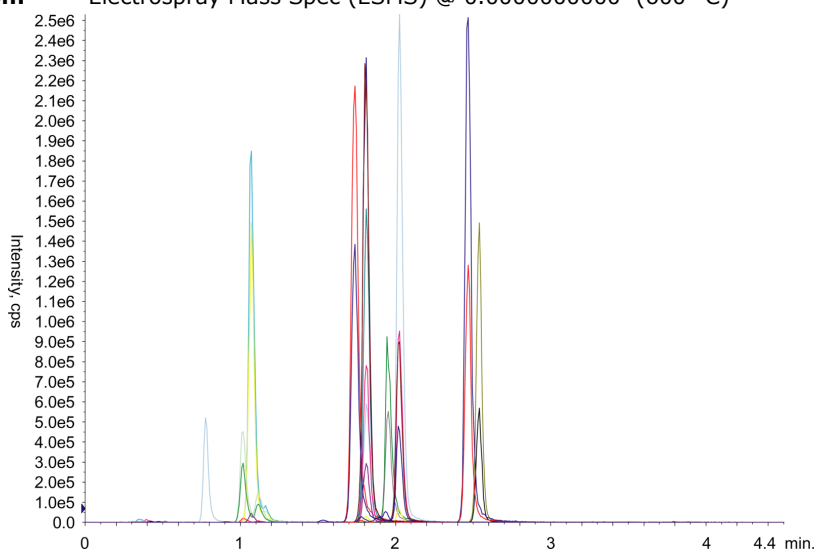


TCA's on Kinetex 2.6u Polar C18, 50x2.1mm (100 ng/mL)**Column:** Kinetex® 2.6 µm Polar C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4759-AN**Elution Type:** Gradient**Eluent A:** 10 mM ammonium formate in water**Eluent B:** 0.1% formic acid in methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	55	45
	2	2	35	65
	3	2.5	5	95
	4	3	5	95
	5	3.1	55	45
	6	4.5	55	45

Flow Rate: 0.5 mL/min**Col. Temp.:** 40 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (600 °C)**ANALYTES:**

- 1 Nortriptyline**
Retention Time: 2.03 min
- 2 Protriptyline**
Retention Time: 1.81 min
- 3 Desmethyldoxepin**
Retention Time: 1.07 min
- 4 Desipramine**
Retention Time: 1.81 min
- 5 Amitriptyline**
Retention Time: 1.95 min
- 6 Doxepin**
Retention Time: 1.01 min
- 7 Imipramine**
Retention Time: 1.74 min
- 8 N-Desmethyldoxepin**
Retention Time: 2.54 min
- 9 Clomipramine**
Retention Time: 2.46 min

©2025 Phenomenex Inc., All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com

support@phenomenex.com